


# Pesticides

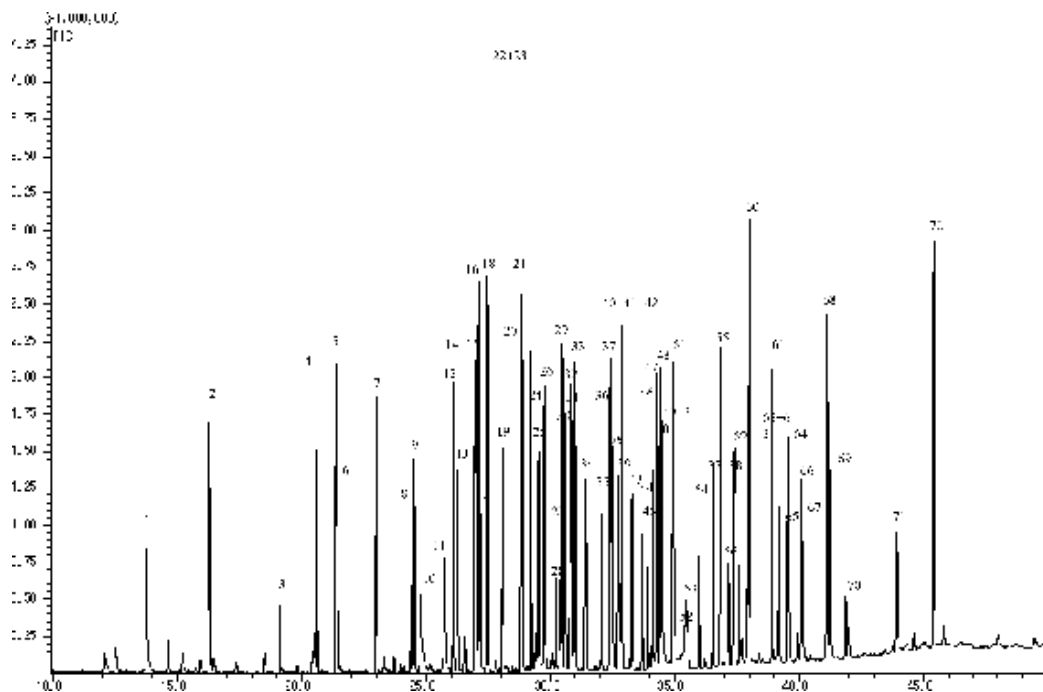

System : GC/MS  
 Column : InertCap 5MS/Sil  
 0.25mm × 30m df=0.25µm  
 Column Temp. : 60°C(1min)→5°C/min→280°C  
 Carrier Gas : He 100kPa  
 Detection : MS

- |           |                      |                    |                        |
|-----------|----------------------|--------------------|------------------------|
| Solutes : | 1. Dichlorvos (DDVP) | 26. Dithiopyr      | 51. Buprofezin         |
|           | 2. Dichlobenil       | 27. Fenitrothion   | 52. Isoxathion         |
|           | 3. Etridiazole       | 28. Probenazole    | 53. Carpropamid        |
|           | 4. Chloroneb         | 29. Esprocarb      | 54. β - Endosulfan     |
|           | 5. Isoprocarb        | 30. Malathion      | 55. Mepronil           |
|           | 6. Molinate          | 31. Chlorpyrifos   | 56. Chlornitrofen      |
|           | 7. Fenobucarb        | 32. Benthocarb     | 57. Propiconazole      |
|           | 8. Trifluralin       | 33. Fenthion       | 58. Edifenphos         |
|           | 9. Benfluralin       | 34. Fthalide       | 59. Endosulfan sulfate |
|           | 10. Pencycuron       | 35. Pendimethalin  | 60. Thenylchlor        |
|           | 11. Dimethoate       | 36. Dimethametryn  | 61. Pyributicarb       |
|           | 12. Simazine         | 37. Isofenphos     | 62. Pyridaphenthion    |
|           | 13. Atrazine         | 38. Methyl dymron  | 63. Iprodione          |
|           | 14. Propyzamide      | 39. Phenthoate     | 64. EPN                |
|           | 15. Pyroquilon       | 40. Procymidone    | 65. Piperophos         |
|           | 16. Diazinon         | 41. Captan         | 66. Anilofos           |
|           | 17. Chlorothalonil   | 42. Dimepiperate   | 67. Bifenox            |
|           | 18. Disulfoton       | 43. Methidathion   | 68. Pyriproxyfen       |
|           | 19. Iprobenfos       | 44. α - Endosulfan | 69. Mefenacet          |
|           | 20. Terbutcarb       | 45. Butamifos      | 70. Benfuracarb        |
|           | 21. Bromobutide      | 46. Napropamide    | 71. Cafenstrole        |
|           | 22. Alachlor         | 47. Flutranil      | 72. Ethofenprox        |
|           | 23. Tolclophosmethyl | 48. Pretilachlor   |                        |
|           | 24. Simetryn         | 49. Isoprothiolane |                        |
|           | 25. Metalaxyl        | 50. Tricyclazole   |                        |