


InertSearch™ for GC

InertCap® Applications

Pesticides

Data No. GA006-0000

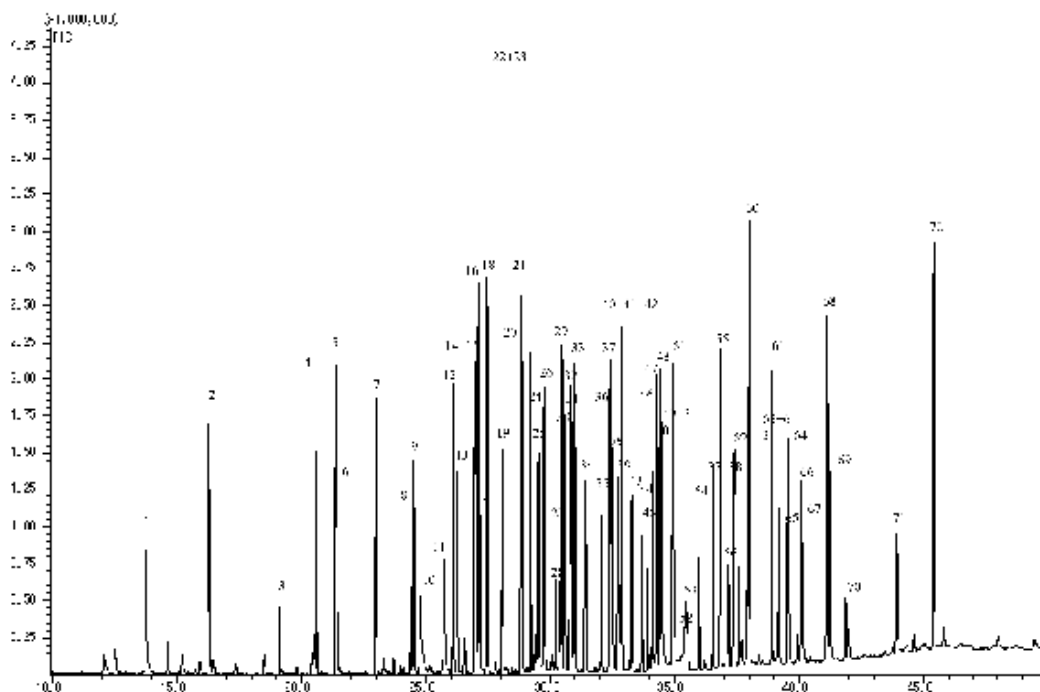

System : GC/MS
Column : InertCap 5MS/Sil
0.25mm×30m d f =0.25µm
Column Temp. : 60°C (1min)→5°C/min→280°C
Carrier Gas : He 100kPa
Detection : MS

- Solutes :
- | | | |
|----------------------|--------------------|------------------------|
| 1. Dichlorvos (DDVP) | 26. Dithiopyr | 51. Buprofezin |
| 2. Dichlobenil | 27. Fenitrothion | 52. Isoxathion |
| 3. Etridiazole | 28. Probenazole | 53. Carpropamid |
| 4. Chloroneb | 29. Esprocarb | 54. β-Endosulfan |
| 5. Isoprocab | 30. Malathion | 55. Mepronil |
| 6. Molinate | 31. Chlorpyrifos | 56. Chlornitrofen |
| 7. Fenobucarb | 32. Benthio carb | 57. Propiconazole |
| 8. Trifluralin | 33. Fenthion | 58. Edifenphos |
| 9. Benfluralin | 34. Fthalide | 59. Endosulfan sulfate |
| 10. Pencycuron | 35. Pendimethalin | 60. Thenylchlor |
| 11. Dimethoate | 36. Dimethametryn | 61. Pyributicarb |
| 12. Simazine | 37. Isofenphos | 62. Pyridaphenthion |
| 13. Atrazine | 38. Methyl dymron | 63. Iprodione |
| 14. Propyzamide | 39. Phenthoate | 64. EPN |
| 15. Pyroquilon | 40. Procymidone | 65. Piperophos |
| 16. Diazinon | 41. Captan | 66. Anilofos |
| 17. Chlorothalonil | 42. Dimepiperate | 67. Bifenox |
| 18. Disulfoton | 43. Methidathion | 68. Pyriproxyfen |
| 19. Iprobenfos | 44. α-Endosulfan | 69. Mefenacet |
| 20. Terbucarb | 45. Butamifos | 70. Benfuracarb |
| 21. Bromobutide | 46. Napropamide | 71. Cafenstrole |
| 22. Alachlor | 47. Flutranil | 72. Ethofenprox |
| 23. Tolclophosmethyl | 48. Pretilachlor | |
| 24. Simetryn | 49. Isoprothiolane | |
| 25. Metalaxyl | 50. Tricyclazole | |