


# Water quality standard restriction pesticides

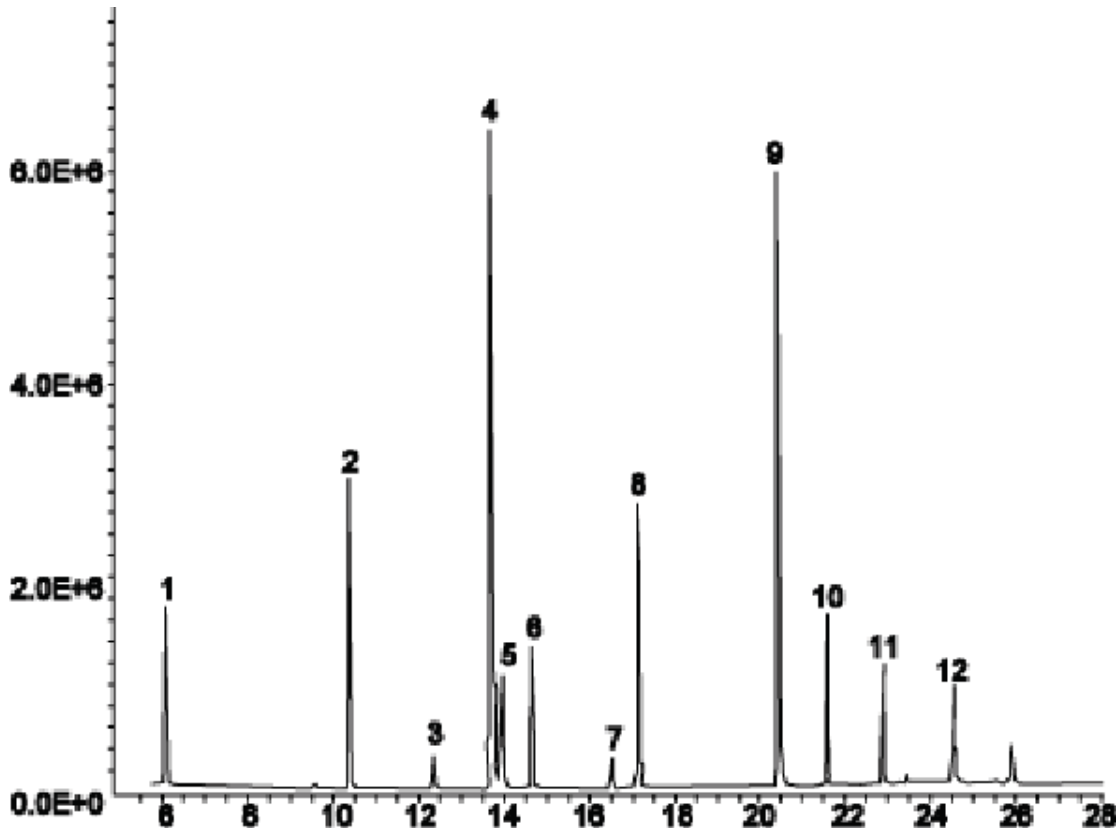

System : GC/MS  
 Column : InertCap 1  
 0.25mmI.D. × 30m df=0.4µm  
 Column Temp. : 50°C(0.5min)→20°C/min→160°C(0.5min)  
 →4°C/min→210°C(6min)→10°C/min→260°C  
 Injection : 1µL  
 Detection : MS  
 Sample : Pesticides

Solutes : 1. Dichlorvos(DDVP)  
 2. Fenobucarb(BPMC)  
 3. Simazine(CAT)  
 4. Chlorothalonil(TPN)+Propyzamide  
 5. Diazinon  
 6. Iprobenfos(IBP)  
 7. Fenitrothion(MEP)  
 8. Thiobencarb  
 9. Isoprothiolane  
 10. Isoxathion  
 11. Chlornitrofen(CNP)  
 12. EPN