


Water quality standard restriction pesticides

Data No. GA007-0000

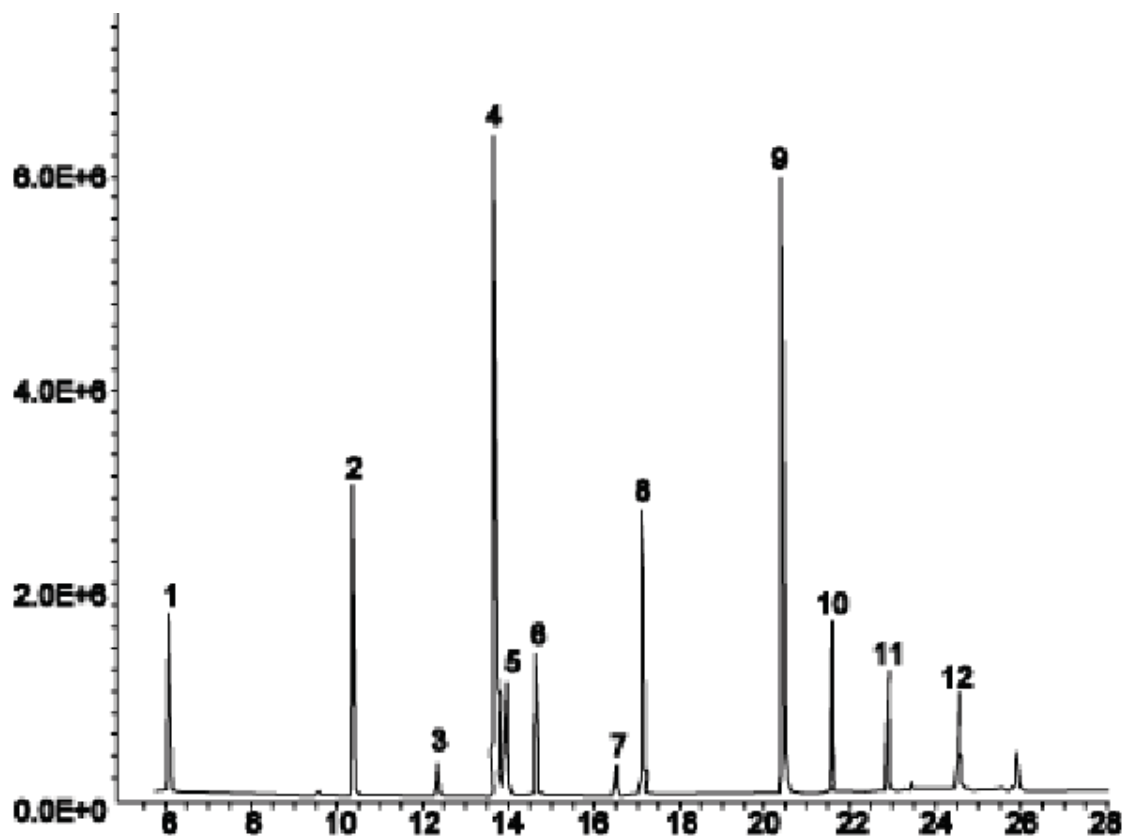

System : GC/MS
Column : InertCap 1
0.25mm I.D. × 30m df=0.4µm
Column Temp. : 50° C (0.5min) → 20° C/min → 160° C (0.5min)
→ 4° C/min → 210° C (6min) → 10° C/min → 260° C
Injection : 1µL
Detection : MS
Sample : Pesticides

Solutes : 1. Dichlorvos(DDVP)
2. Fenobucarb(BPMC)
3. Simazine(CAT)
4. Chlorothalonil(TPN)+ Propyzamide
5. Diazinon
6. Iprobenfos(IBP)
7. Fenitrothion(MEP)
8. Thiobencarb
9. Isoprothiolane
10. Isoxathion
11. Chlornitrofen(CNP)
12. EPN