


# Pesticides in golf courses

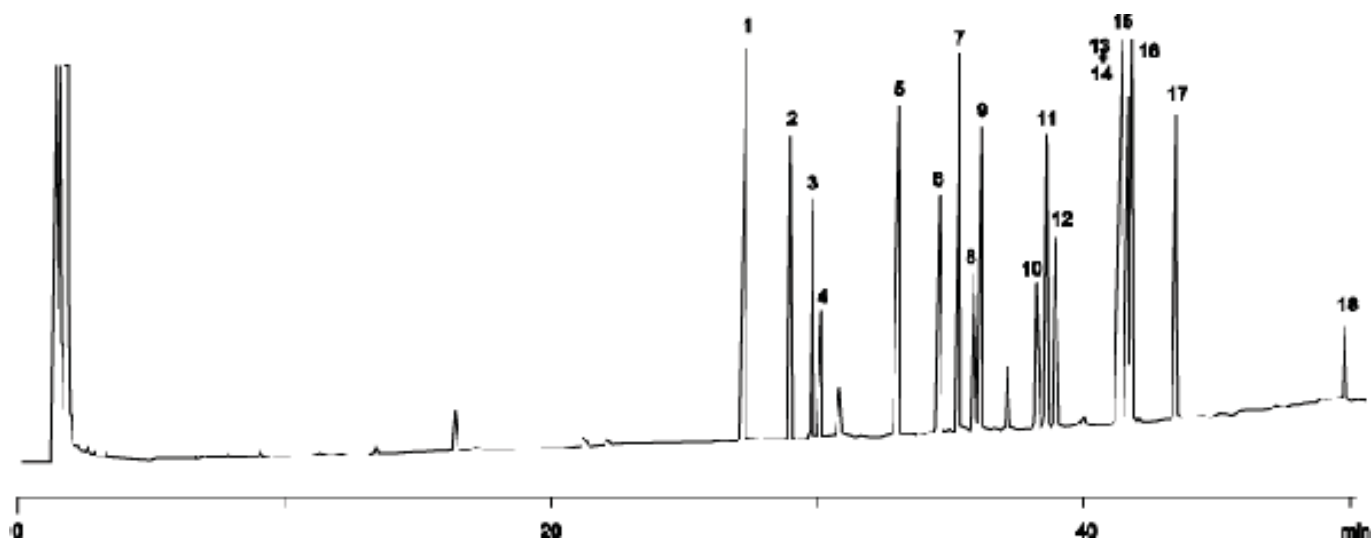

System : GC/FID  
 Column : InertCap 5  
           0.25mmI.D. × 30m  
 Column Temp. : df=0.4μm  
 Carrier Gas : 100°C → 3°C/min → 250°C  
 Injection : He 30cm/sec  
 Detection : Split 1 : 70  
 Sample : FID Range 1 × 8  
           : Pesticides

Solutes : 1. Simazine(CAT)  
           2. Propyzamide  
           3. Diazinon  
           4. Chlorothalonil(TPN)  
           5. Tolclophosmethyl  
           6. Fenitrothion(MEP)  
           7. Thiobencarb  
           8. Fenthion(MPP)  
           9. Chlorpyrifos  
           10. Pendimethalin  
           11. Captan  
           12. Isofenphos  
           13. Napropamide  
           14. Butamifos  
           15. Flutoluanil  
           16. Isoprothiolane  
           17. Isoxathion  
           18. Iprodione