


## Organophosphorus pesticides

Data No. GA010-0000

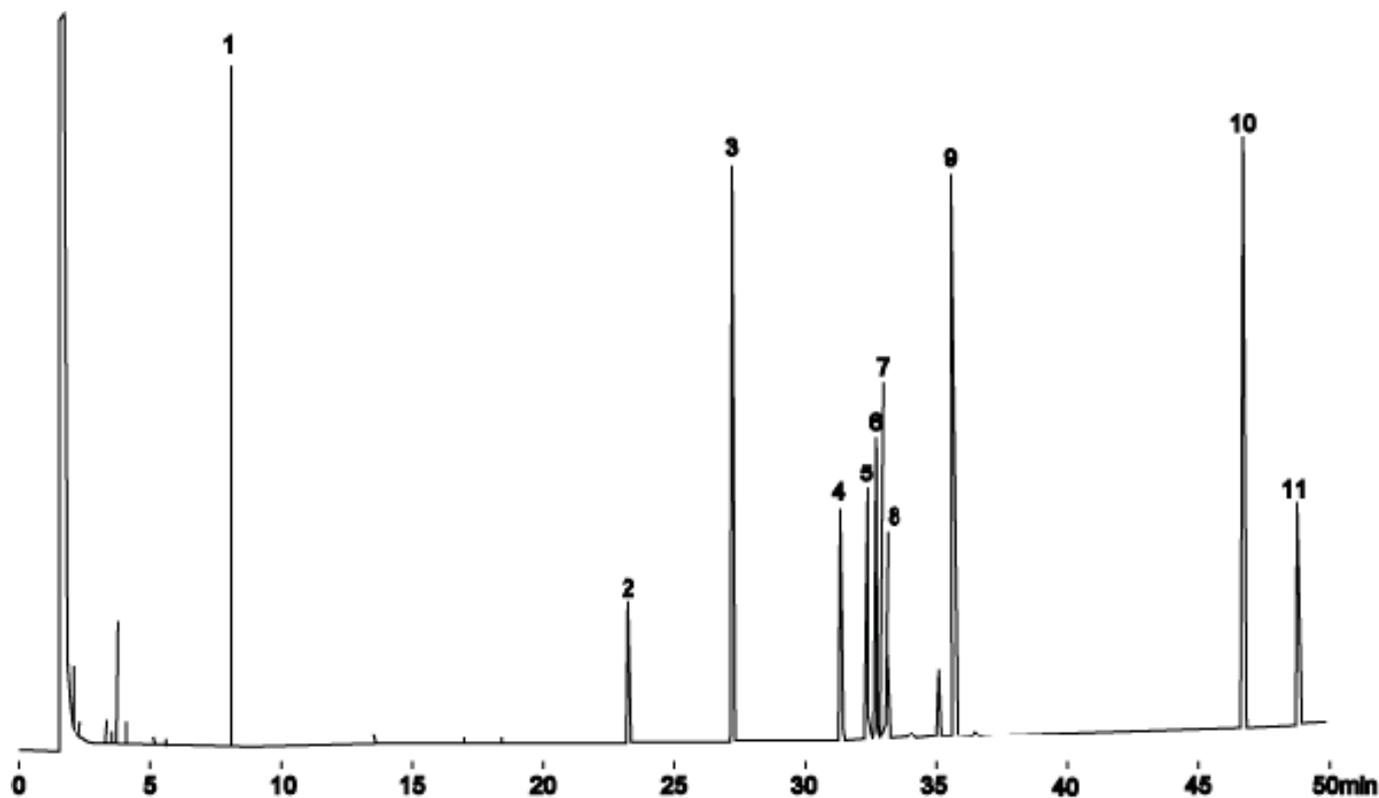

System : GC/FID  
Column : InertCap 1  
0.25mmI.D.×30m df=0.4µm  
Column Temp. : 100°C→3°C/min→270°C  
Carrier Gas : He 80kPa  
Injection : Split 1 : 50 1µl  
Detection : FID Range 1×256  
Sample : Pesticides

Solutes : 1. Dichlorovos  
2. Dimethoate  
3. Diazinon  
4. Fenitrothion  
5. Malathion  
6. Fenthion  
7. Parathion  
8. Chlorpyriphos  
9. Chlorfenvinphos (CVP)  
10. Phenthoate (PAP)  
11. EPN  
12. Phosalone