


Indoor air by AUTOCAN

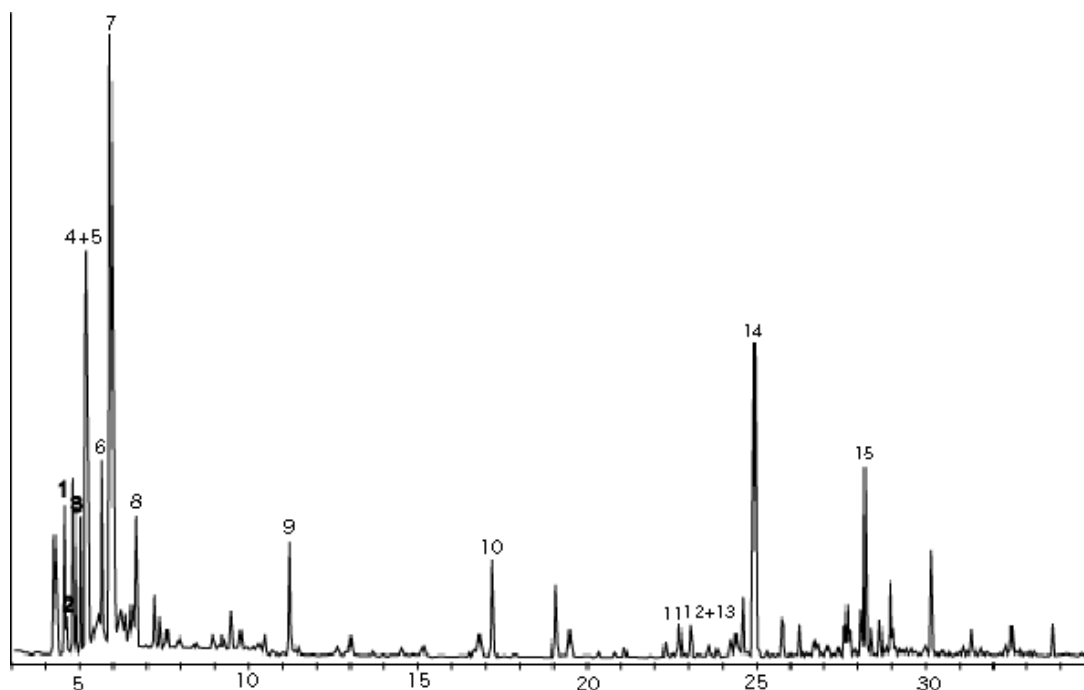

System : GC/MS AUTO Can
 Column : InertCap AQUATIC
 0.25mmI.D. × 60m df=1.0µm
 Carrier Gas : He
 Injection : AUTO Can/Cryofocus 200mL
 Detection : MS
 Sample : Indoor air

Solutes : 1. Propane
 2. Unknown
 3. Butane
 4. Methanol
 5. Acetaldehyde
 6. 2-Methylbutane
 7. Ethanol
 8. Acetone
 9. Ethylacetate
 10. Toluene
 11. Ethylbenzene
 12. m-Xylene
 13. p-Xylene
 14. α-Pinene
 15. Limonene