


# InertSearch™ for GC

InertCap® Applications

## VOC & Amines

Data No. GA065-0000

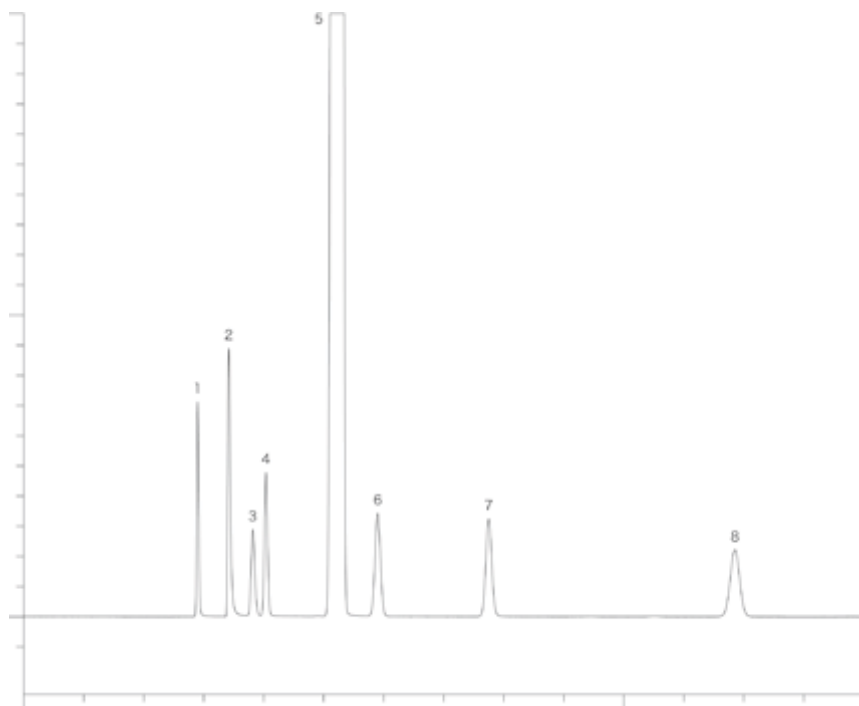

System : GC/FID  
Column : InertCap for Amines  
0.32mmI.D. × 30m  
Column Temp. : 50°C  
Carrier Gas : He 70kPa  
Injection : Split 1 : 50  
Detection : FID

Solutes : 1. Methanol  
2. Dimethylamine (DMA)  
3. Trimethylamine (TMA)  
4. Ethanol  
5. iso-Propanol  
6. n-Pentane  
7. n-Propanol  
8. n-Hexane