


# InertSearch™ for GC

InertCap® Applications

## Cresols

Data No. GA070-0000

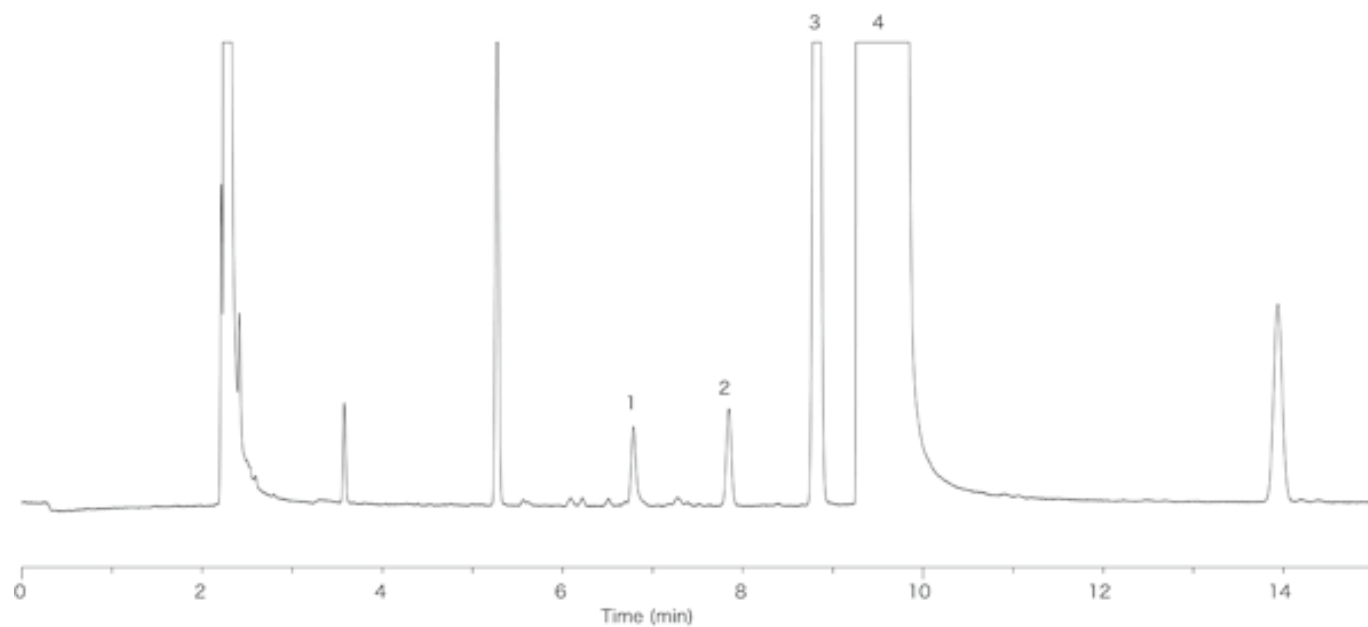

System : GC/FID  
Column : InertCap CHIRAMIX  
0.25mmI.D.×30m df=0.25µm  
Column Temp. : 150°C  
Carrier Gas : He 100kPa  
Injection : Split 1 : 50 0.5µL  
250°C  
Detection : FID  
250°C  
Sample : Cresol

Solutes : 1. Phenol  
2. o-Cresol  
3. p-Cresol  
4. m-Cresol