


InertSearch™ for GC

InertCap® Applications

Glycerin, Diethylene glycol

Data No. GA071-0000

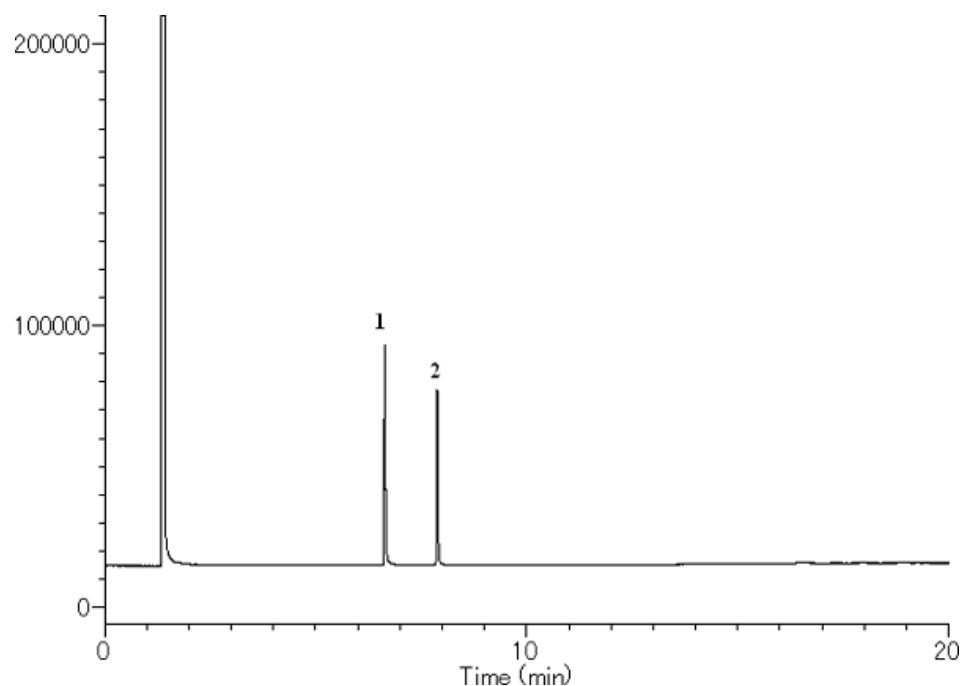

System : GC/FID
Column : InertCap1701 0.32mmI.D. ×30m df=1.0µm
Column Temp. : 100°C→7.5°C→220°C(4min hold)
Carrier Gas : 100kPa
Injection : Split1:20 220°C 1µL
Detection : FID Range10⁰ 250°C
Sample : each 0.5mg/mL

Solutes : 1. Diethylene glycol
2. Glycerin