


InertSearch™ for GC

InertCap® Applications

Isomers of cresol

Data No. GA083-0000

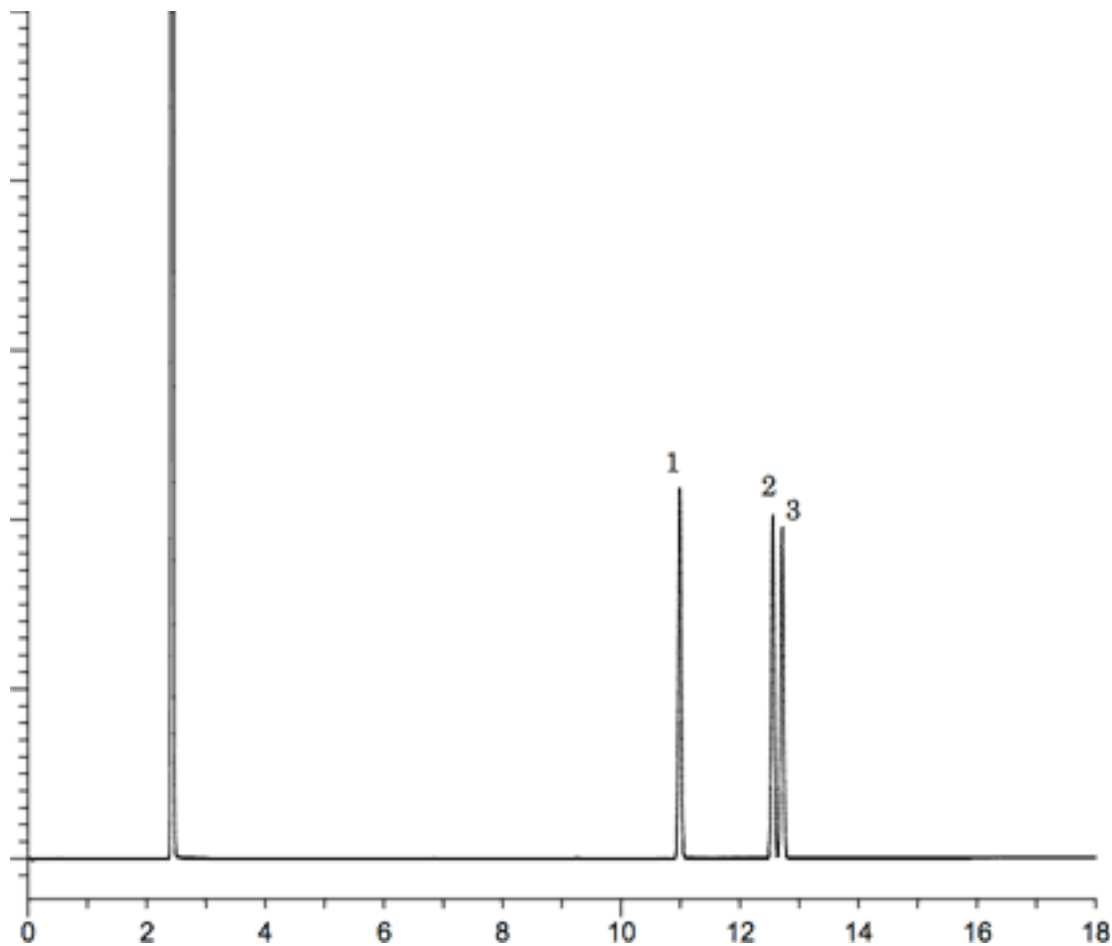

System : GC/FID
Column : InertCap Pure-WAX
0.25mm I.D. × 30m df=0.25µm
Column Temp. : 150°C (5min) → 5°C/min → 250°C
Carrier Gas : He 100kPa
Injection : Split 1 : 50 0.2µl
250°C
Detection : FID Range 10°
260°C
Sample : Cresol isomers

Solutes : 1. o-Cresol
2. p-Cresol
3. m-Cresol