


## Residual solvents in pharmaceuticals Class1

Data No. GA090-0000

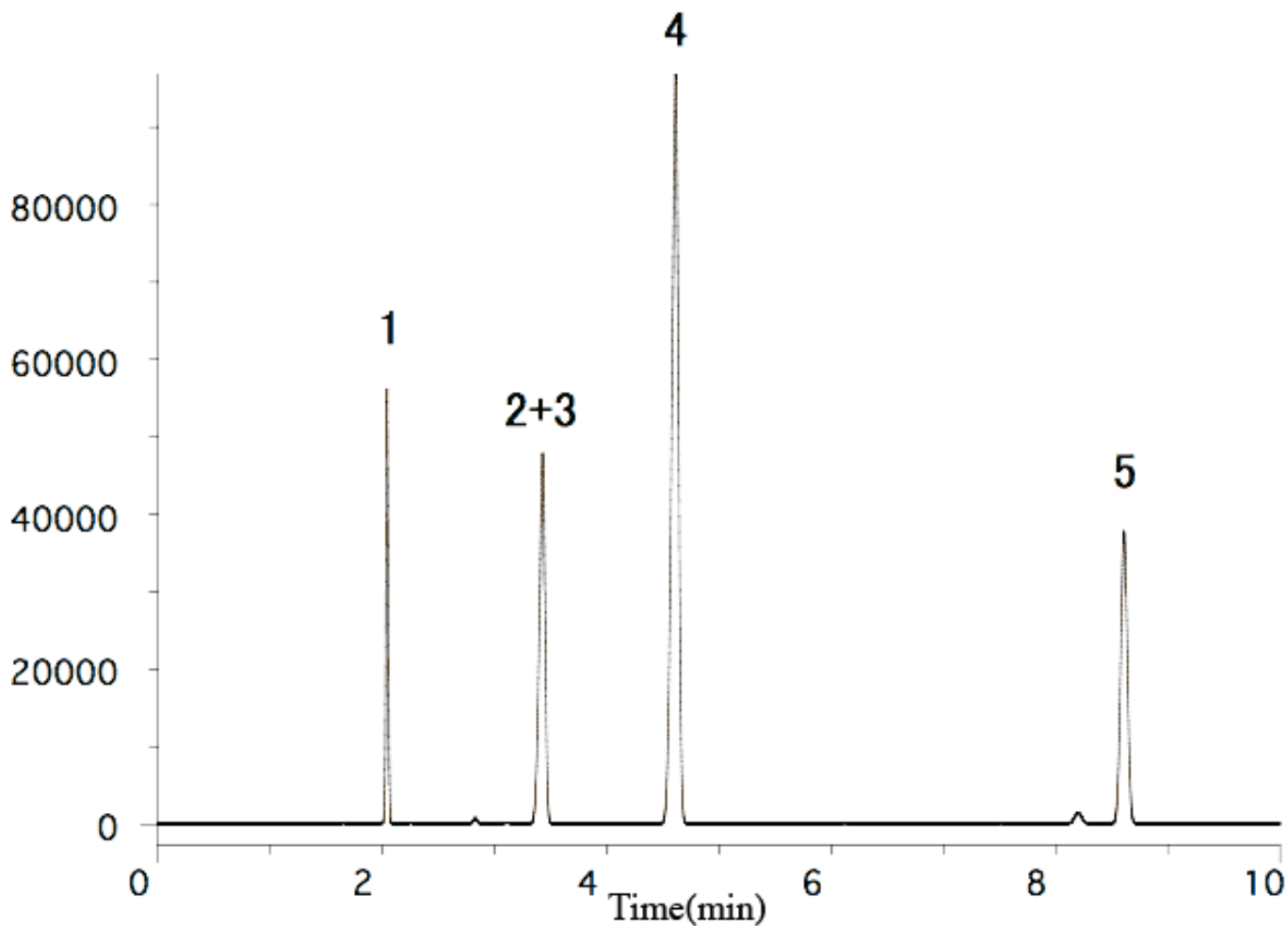

System : GC/FID  
Column : InertCap Pure-WAX  
0.53mm I.D. × 30m df=1µm  
Column Temp. : 40°C – 5°C/min – 90°C  
Carrier Gas : He 20kPa 35cm/sec  
Injection : Split Flow 200ml/min  
240°C  
Detection : FID Range 10<sup>2</sup>  
240°C  
Sample : Solvents  
Mixed evenly  
Sample Size : 0.4µL

Solutes 1. 1,1-Dichloroethylene  
2. 1,1,1-Trichloroethane  
3. Carbon tetrachloride  
4. Benzene  
5. 1,2-Dichloroethane