


Glycols & Glycerine

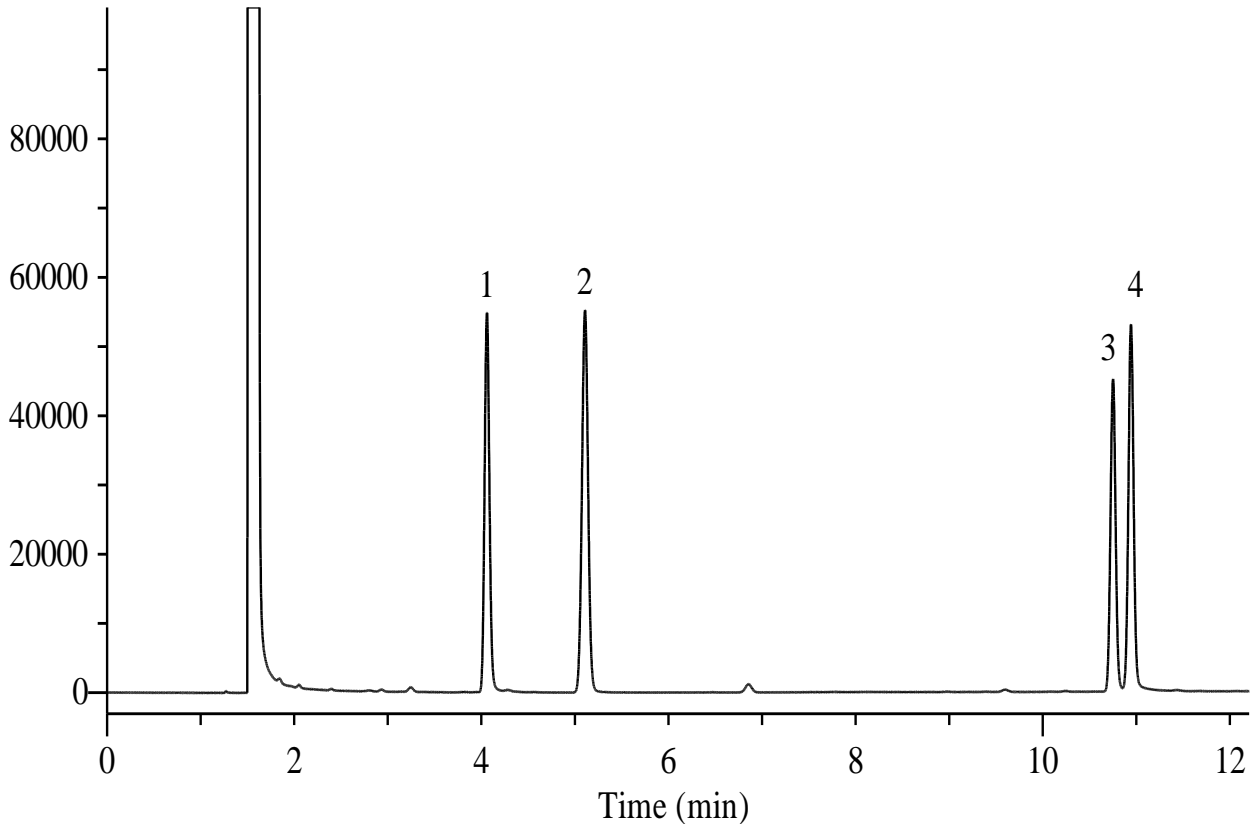


Conditions

System : GC/FID
Column : InertCap 1
 0.32 mm I.D. x 30 m df = 5.0 µm
Col. Cat. No. : 1010-11249
Col. Temp. : 100 °C (5 min hold) - 10 °C/min - 240 °C
Carrier Gas : He 100 kPa
Injection : Split flow 53.6 mL/min
 240 °C
Detection : FID Range 10⁰
 240 °C
Sample Size : 500 µg/mL in Methanol
 1 µL


Analyte : 1. Ethylene glycol
 2. Propylene glycol
 3. Glycerine
 4. Diethylene glycol

Glycols & Glycerine


Chemical Structures

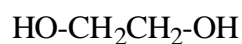


1. Ethylene glycol

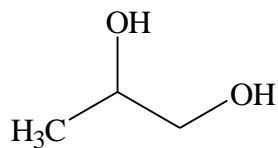

2. Propylene glycol

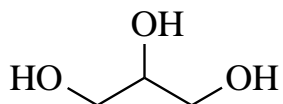

3. Glycerine

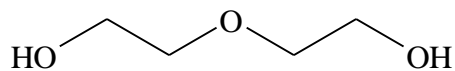

4. Diethylene glycol