


# InertSearch™ for GC

InertCap® Applications

## Organochlorine Pesticide

Data No. GA139-0875

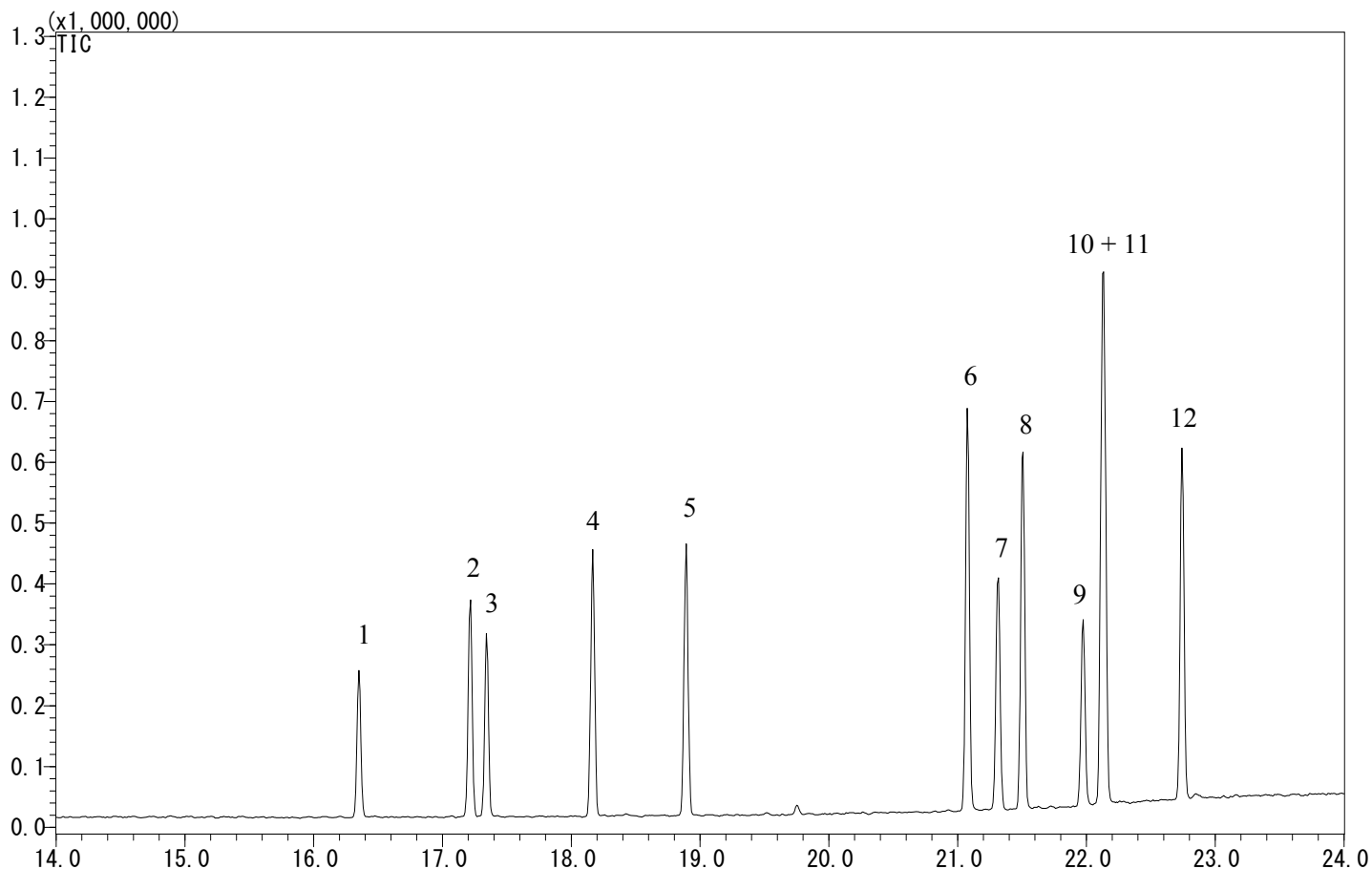

### Conditions

**System** : GC/MS  
**Column** : InertCap 35  
0.25 mm I.D. x 30 m df = 0.25 µm  
**Col. Cat. No.** : 1010-63142  
**Col. Temp.** : 60 °C - 10 °C/min - 290 °C (7 min hold)  
**Carrier Gas** : He 35 cm/sec  
**Injection** : Split 1:30  
250 °C  
**Detection** : MS Scan (45 - 500m/z)  
Interface Temp. 280 °C  
**Sample Size** : 100 µg/mL in Isooctane  
1 µL

**Analyte** : 1. α-HCH  
2. γ-HCH  
3. β-HCH  
4. Heptachlor  
5. Aldrin  
6. p,p'-DDE  
7. Dieldrin  
8. o,p'-DDD  
9. Endrin  
10. p,p'-DDD  
11. o,p'-DDT  
12. p,p'-DDT