
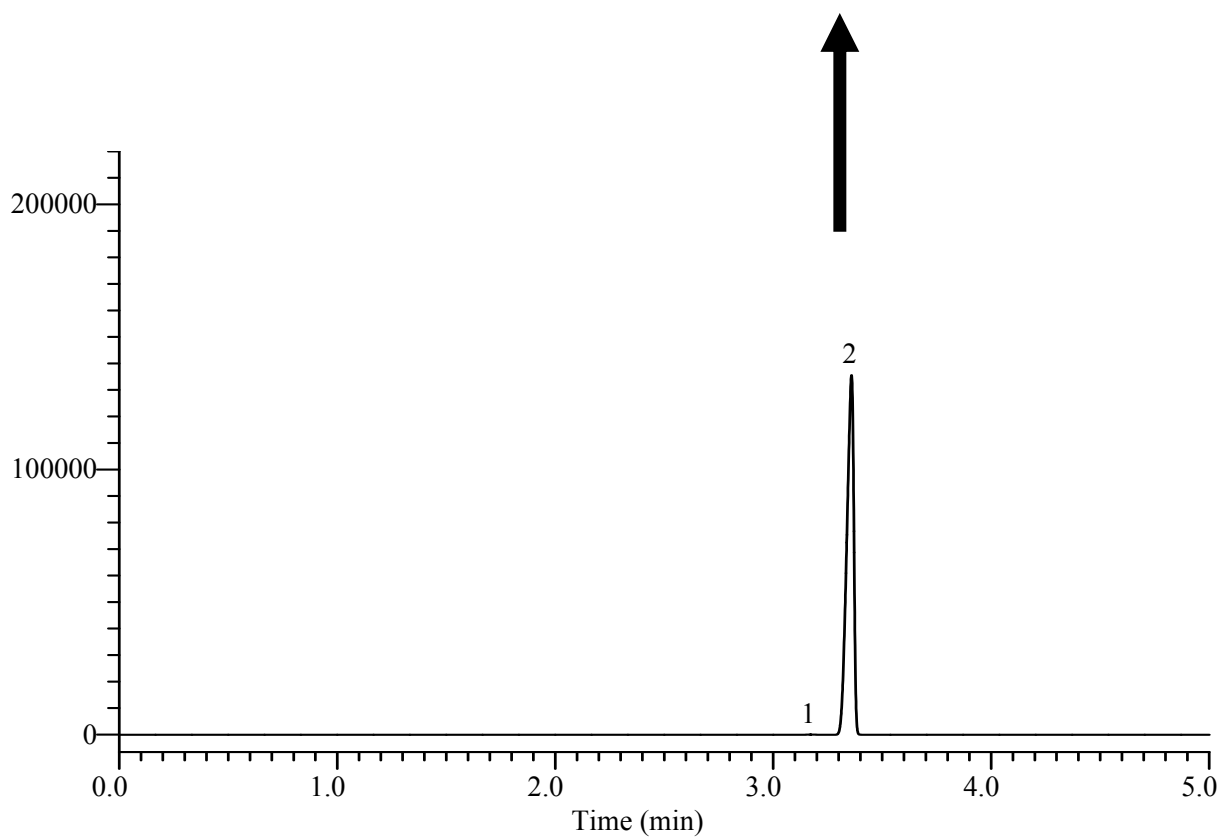


# InertSearch™ for GC

InertCap® Applications

## Impurities in 2-Propanol (Isopropyl alcohol)

Data No. GA159-0591

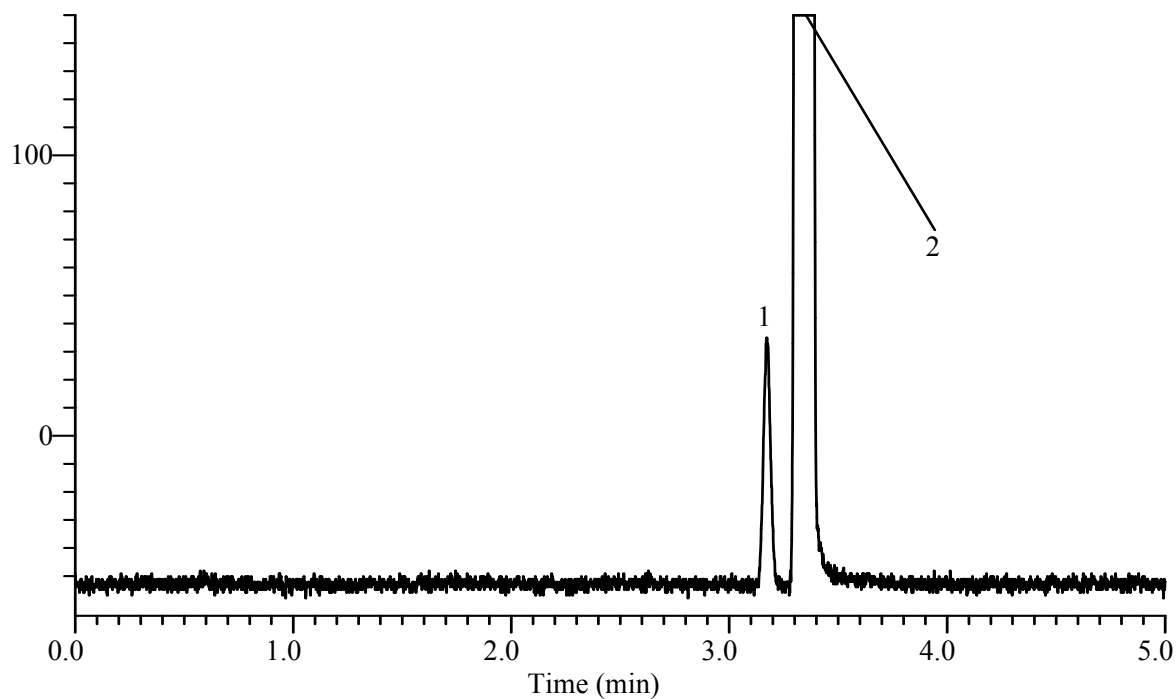

# InertSearch™ for GC

InertCap® Applications

## Impurities in 2-Propanol (Isopropyl alcohol)

Data No. GA159-0591

### Conditions

**System** : GC/FID  
**Column** : InertCap WAX  
0.25 mm I.D. x 30 m df = 0.25 µm  
**Col. Cat. No.** : 1010-67142  
**Col. Temp.** : 70 °C  
**Carrier Gas** : He 70 kPa  
**Injection** : Split flow 180 mL/min  
150 °C  
**Detection** : FID Range 10<sup>2</sup>  
150 °C  
**Sample Size** : 1 mg/mL in 2-Propanol  
0.5 µL

### Analyte

1. Methanol
2. 2-Propanol (Isopropyl alcohol)