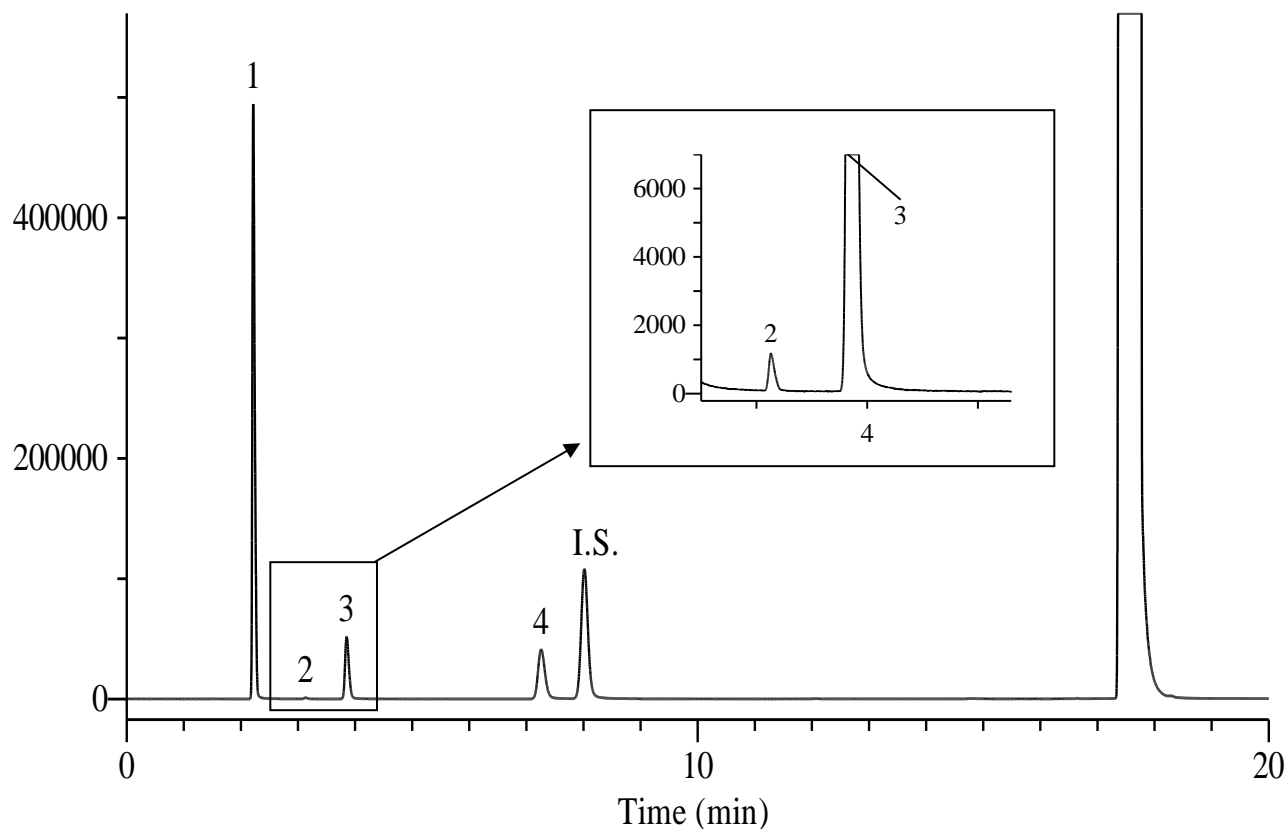


Purity test of Epirubicin Hydrochloride by Japanese Pharmacopoeia


Conditions

System : GC/FID
Column : InertCap Pure-WAX
 0.53 mm I.D. x 30 m df = 1.00 μm
Col.Cat. No. : 1010-68445
Col. Temp. : 40 °C (11 min hold) - 10 °C/min
 - 90 °C - 50 °C/min - 130 °C (30 min hold)
Carrier Gas : He 20 kPa
Injection : Split 1:15
 120 °C
Detection : FID Range 10¹⁰
 150 °C
Sample Size : Analyte in *N,N*-Dimethylformamide
 1.0 μL

Analyte : 1. Acetone 1.25 μL/mL
 2. Methanol 0.0068 μL/mL
 3. Ethanol 0.30 μL/mL
 4. 1-Propanol 0.32 μL/mL
 I.S. 1,4-Dioxane 1.00 μL/mL
 *Standard mixture
 R(1-I.S.) : 41.0 > 30(JP15)

Chemical Structure
