


InertSearch™ for GC

InertCap® Applications

ε-Caprolactam

Data No. GA171-0904

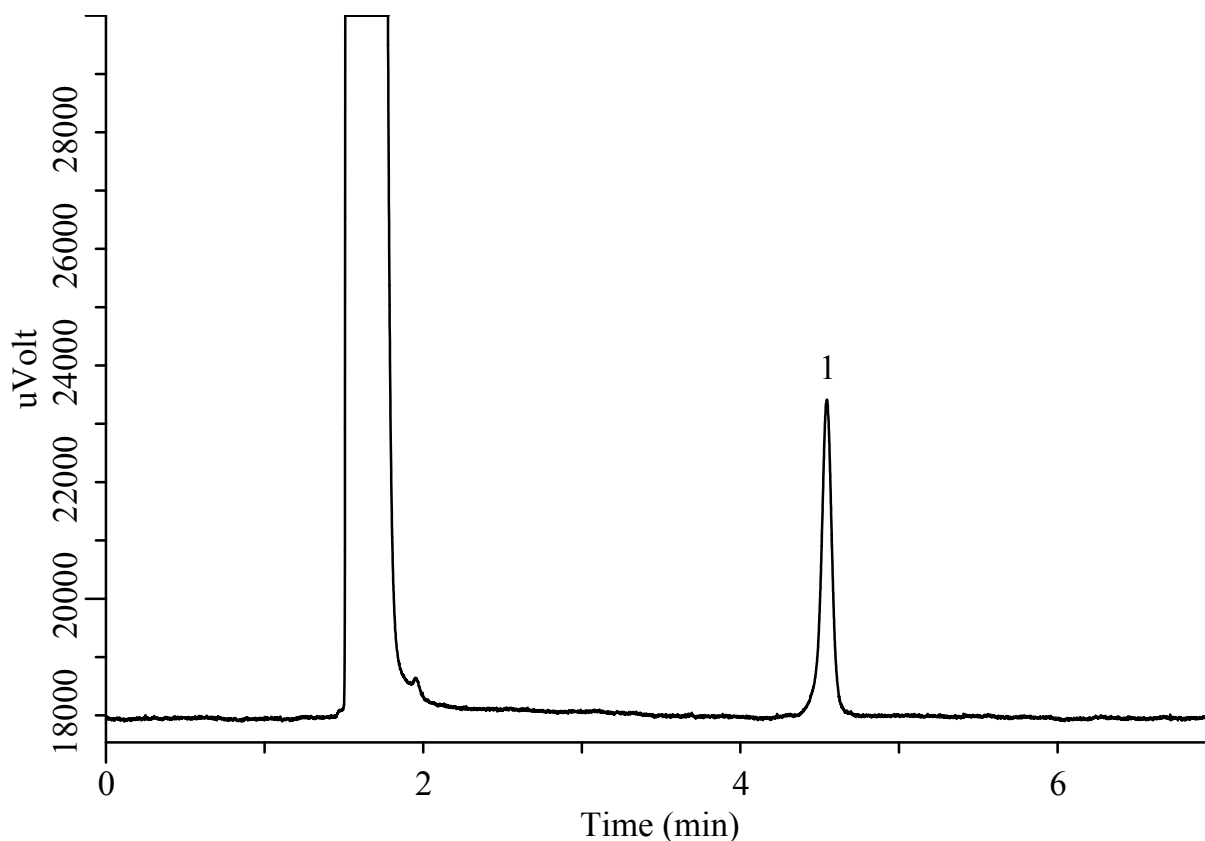

Conditions

System : GC/FID
Column : InertCap 1
0.32 mm I.D. x 30 m df = 5.00 µm
Col.Cat. No. : 1010-11249
Col. Temp. : 240 °C
Carrier Gas : He 110 kPa
Injection : Split flow 18 mL/min
240 °C
Detection : FID Range 10⁰
240 °C
Sample Size : 15 mg/L in 20% Ethanol solution
1.0 µL
Analyte : 1. ε-Caprolactam