


Acrylonitrile

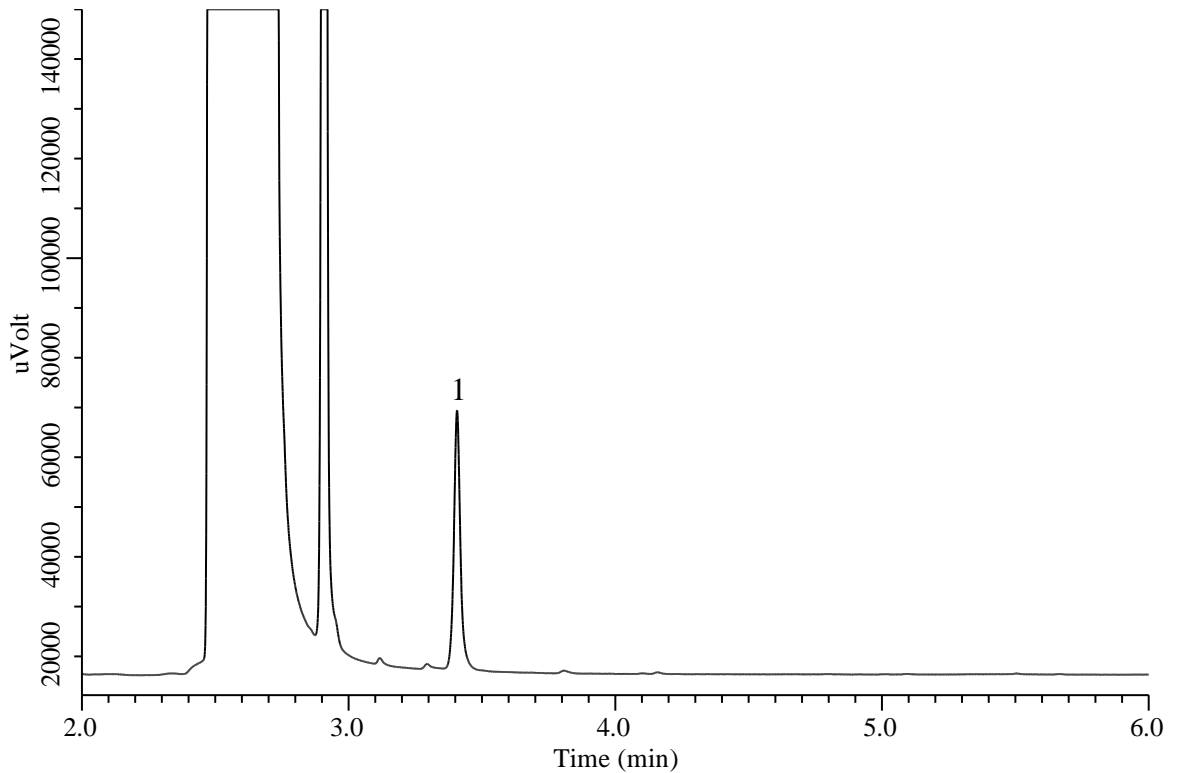

Conditions

System	: GC/FID	Analyte	:
Column	: InertCap Pure-WAX 0.25 mm I.D. x 30 m df = 0.25 µm	1. Acrylonitrile	
Col.Cat. No.	: 1010-68142		
Col. Temp.	: 50 °C (1 min hold) – 10 °C/min – 200 °C		
Carrier Gas	: He 100 kPa		
Injection	: Split flow 12.7 mL/min 200 °C		
Detection	: FID Range 10 ⁰ 250 °C		
Sample Size	: 50 µg/mL in 20% Ethanol solution 1.0 µL		