


Hydrocarbons C1-C2

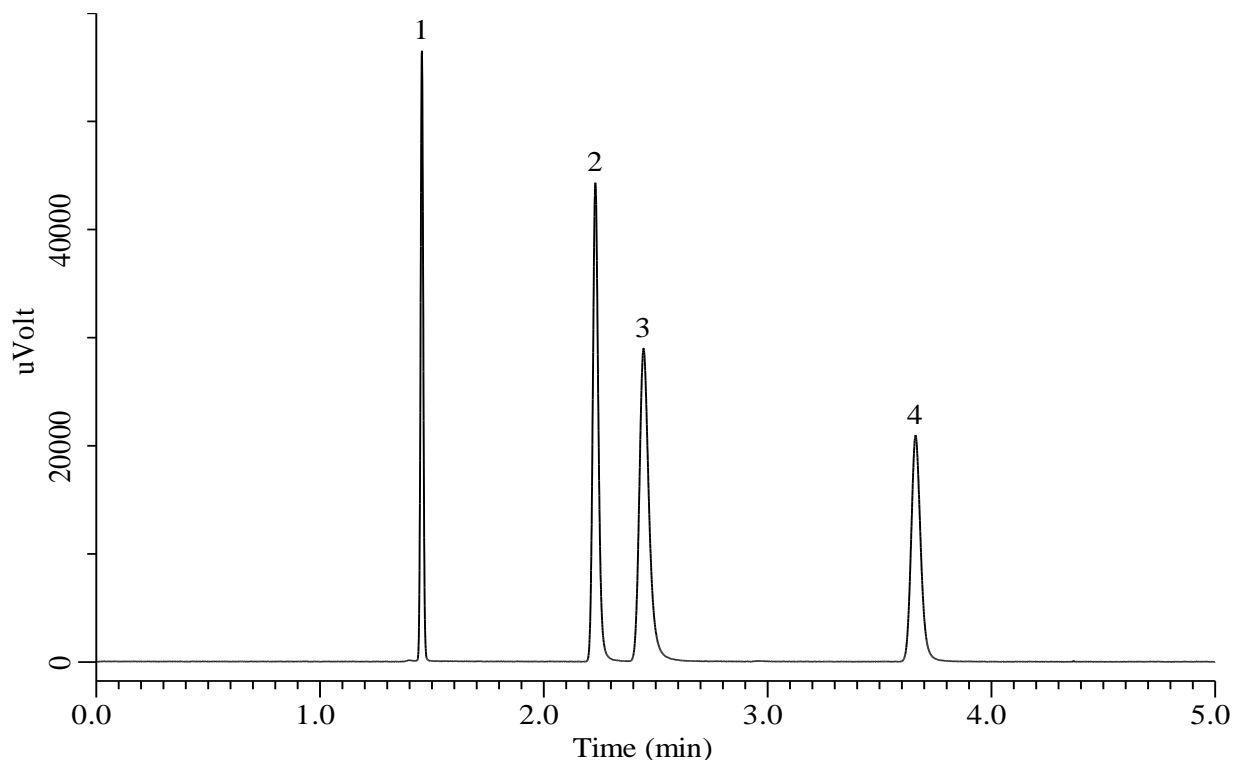

Conditions

System : GC/FID
Column : TC-BOND U
 0.53 mm I.D. x 30 m df = 20 µm
Col. Cat. No. : 1010-40443
Col. Temp. : 40 °C
Carrier Gas : He 30 kPa
Injection : Split flow 50 mL/min
 190 °C
Detection : FID Range 10⁰
 190 °C
Sample Size : 100 ppm Vol/Vol in He
 200 µL

Analyte

1. Methane
 2. Ethylene
 3. Ethane
 4. Acetylene