

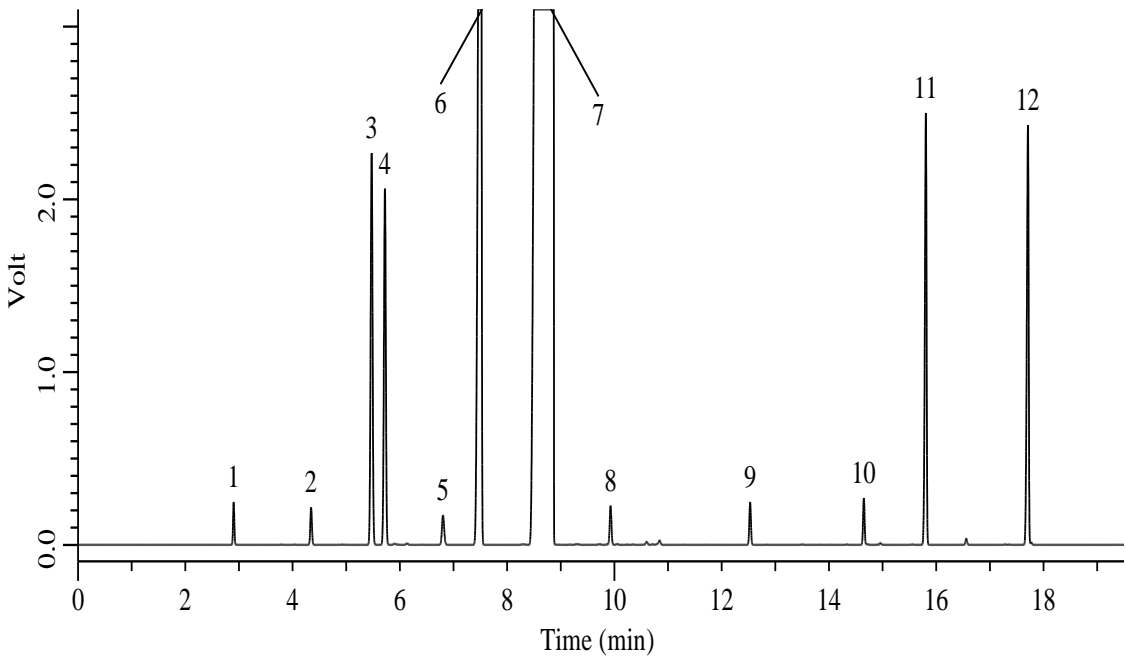
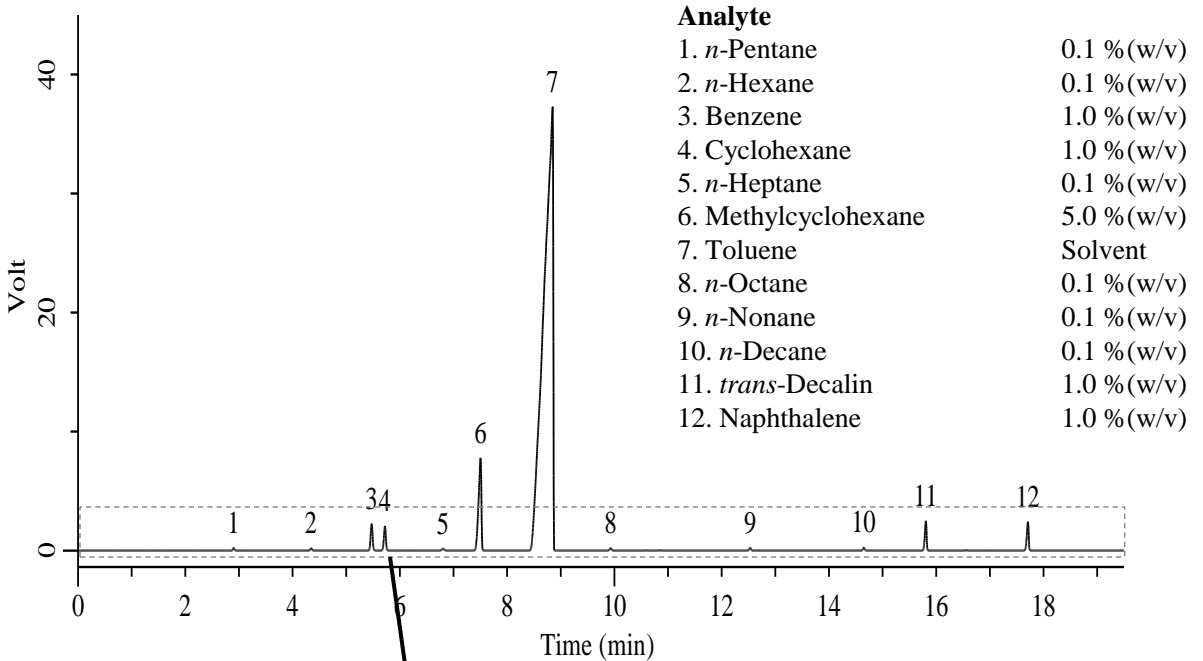


Organic hydride

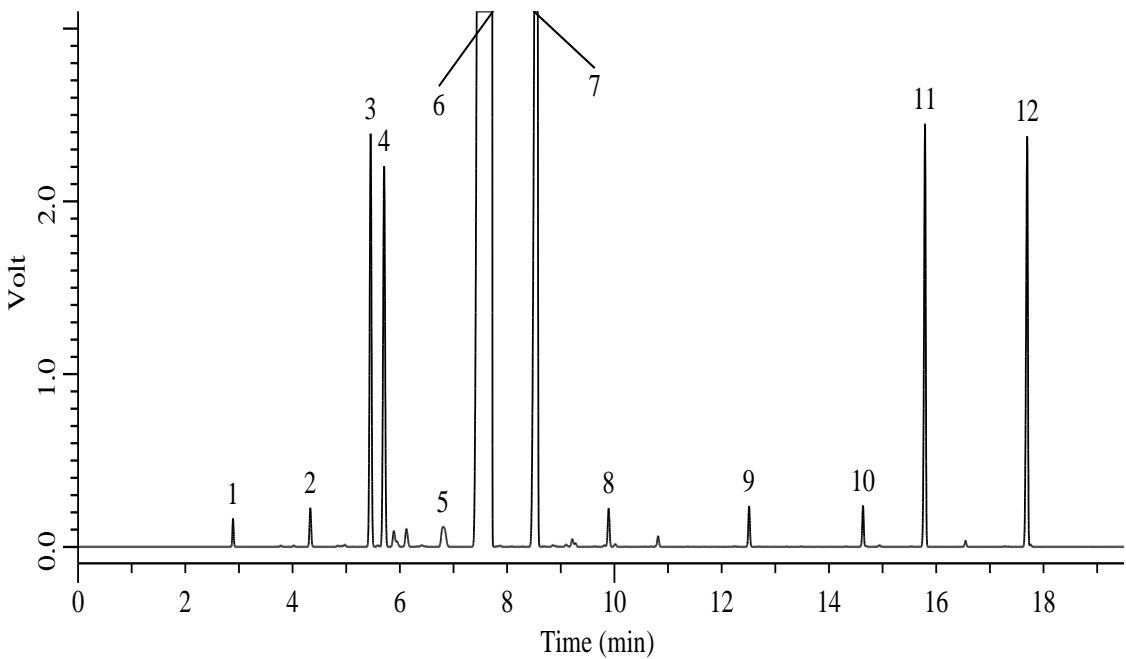

Organic hydride

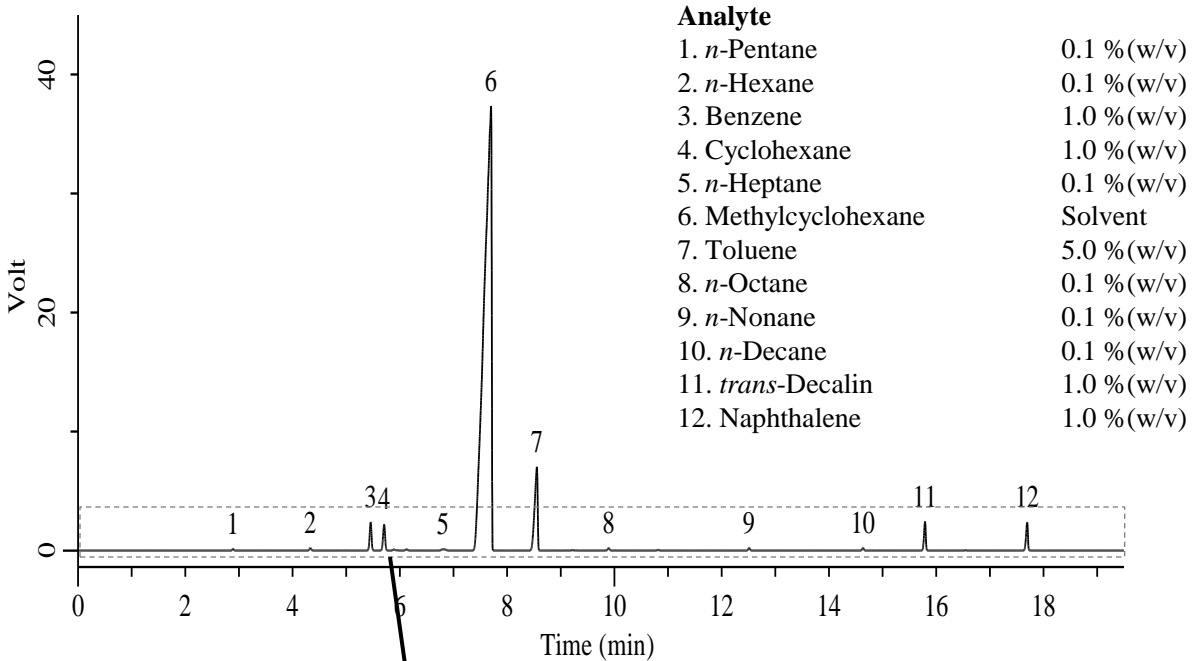

Organic hydride

Conditions

System	: GC-FID
Column	: InertCap 1 0.32 mm I.D. x 30 m df = 1.0 μ m
Col. Cat. No.	: 1010-11245
Col. Temp.	: 40 °C (2 min hold) – 5 °C/min – 75 °C – 10 °C/min – 180 °C
Carrier Gas	: He 1.7 mL/min
Injection	: Split flow 34 mL/min 230 °C
Detection	: FID Auto Range 250 °C
Sample Size	: 0.5 μ L