


Residual solvents in pharmaceuticals (Class 2A)

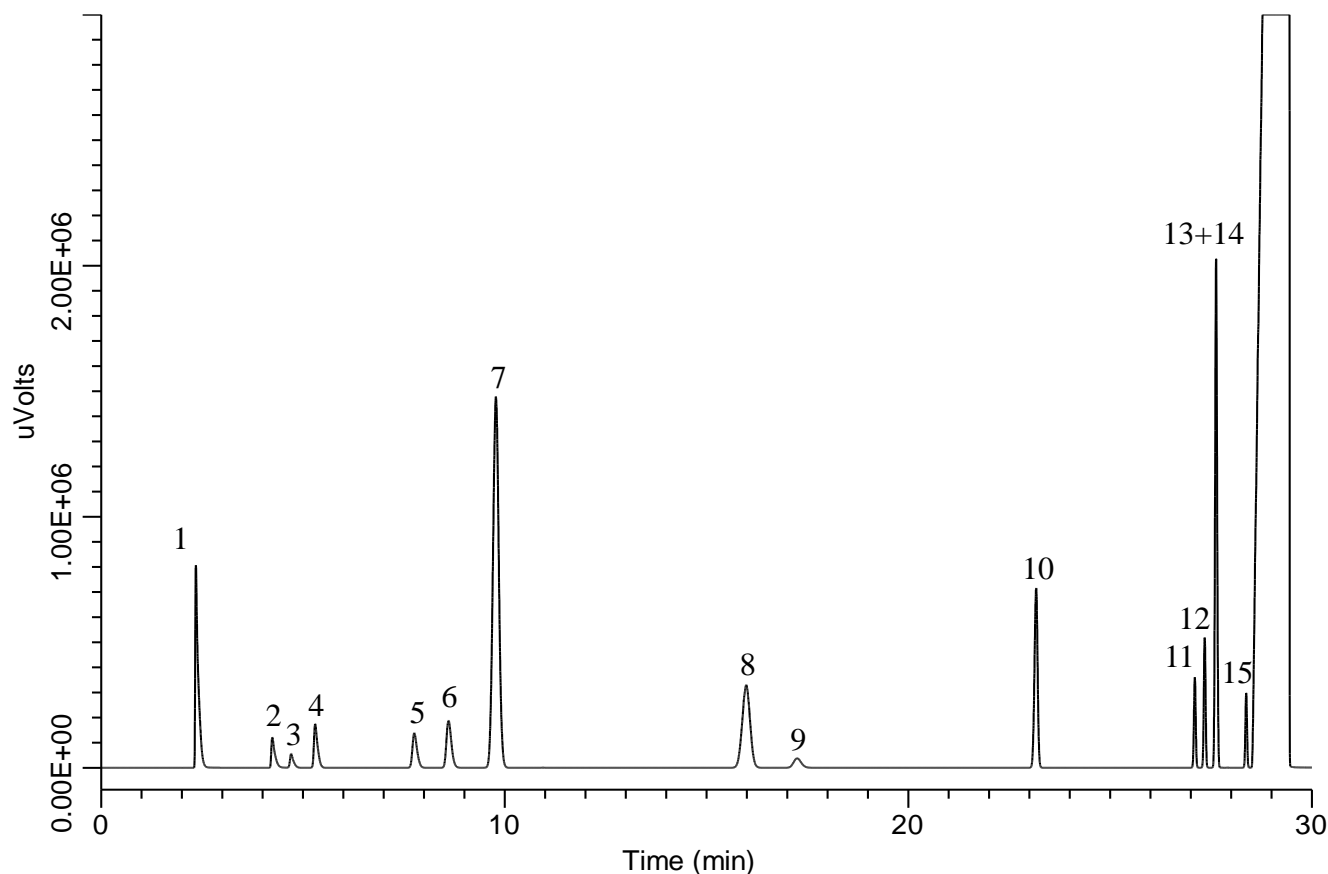

Conditions

System	: GC4000 Plus-FID
Column	: InertCap 624MS 0.32 mm I.D. × 30 m df = 1.80 μm
Col. Cat. No.	: 1010-64747
Col. Temp.	: 40 °C(20 min) - 10 °C/min - 240 °C(20 min)
Carrier Gas	: He 2.20 mL/min
Injection	: Split 20:1 140 °C
Detection	: FID Auto Range 250 °C
Sample Size	: 1.0 μL Analyte in Dimethyl sulfoxide

Analyte

1. Methanol	(15.0 mg/mL)
2. Acetonitrile	(2.05 mg/mL)
3. Dichloromethane	(3.00 mg/mL)
4. <i>trans</i> -1,2-Dichloroethylene	(4.70 mg/mL)
5. <i>cis</i> -1,2-Dichloroethylene	(4.70 mg/mL)
6. Tetrahydrofuran	(3.45 mg/mL)
7. Cyclohexane	(19.4 mg/mL)
8. Methylcyclohexane	(5.90 mg/mL)
9. 1,4-Dioxane	(1.90 mg/mL)
10. Toluene	(4.45 mg/mL)
11. Chlorobenzene	(1.80 mg/mL)
12. Ethylbenzene	(1.84 mg/mL)
13. <i>m</i> -Xylene	(6.51 mg/mL)
14. <i>p</i> - Xylene	(1.52 mg/mL)
15. <i>o</i> -Xylene	(0.98 mg/mL)