


## Flavor by Japanese Standards of Food Additives

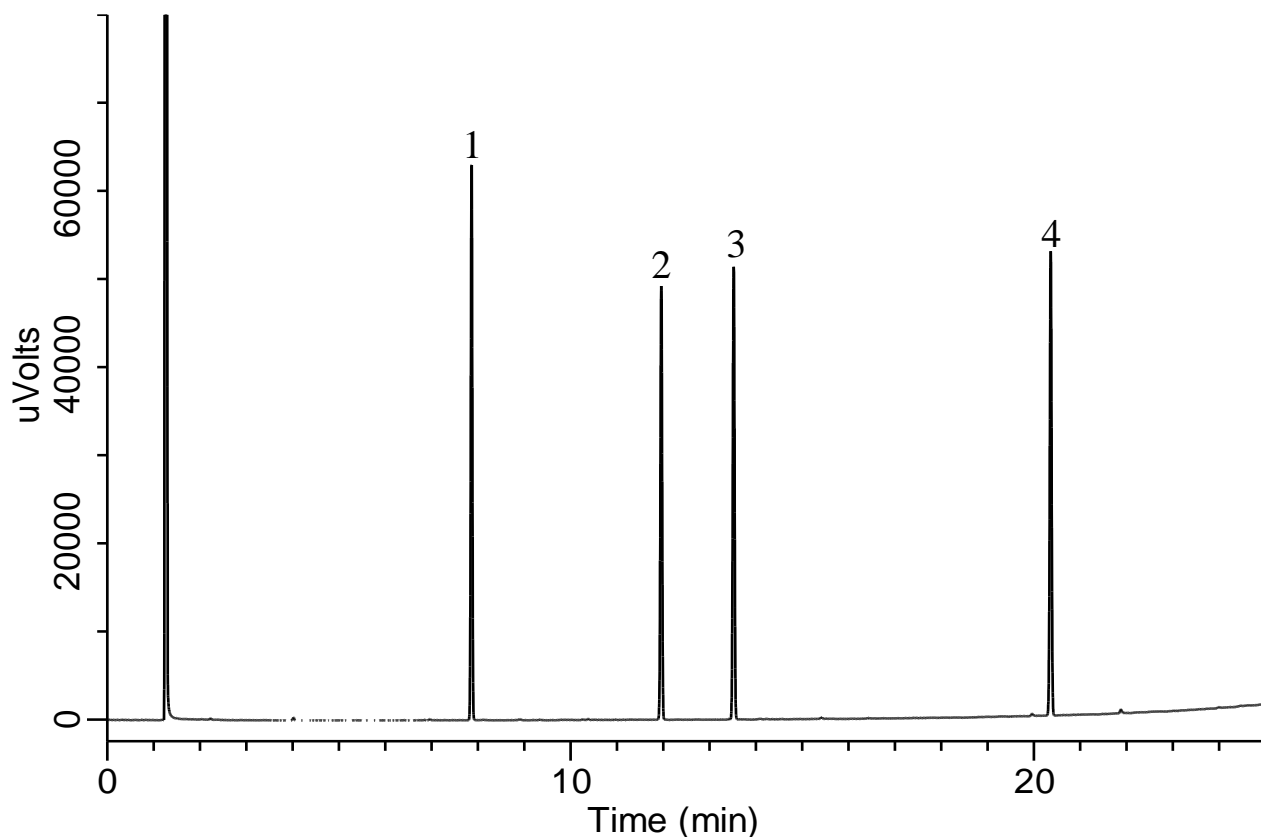

### Conditions

**System** : GC-FID  
**Column** : InertCap Pure-WAX  
           0.25 mm I.D. × 30 m df = 0.25 μm  
**Col. Cat. No.** : 1010-68142  
**Col. Temp.** : 100 °C - 5 °C/min - 230 °C  
**Carrier Gas** : He 200 kPa  
**Injection** : Split flow 320 mL/min  
           250 °C  
**Detection** : FID Auto Range  
           250 °C  
**Sample Size** : 1000 μg/mL in Methanol  
           1 μL

### Analyte

1. 5,6,7,8-Tetrahydroquinoxaline  
 2. 5-Methylquinoxaline  
 3. 6-Methylquinoline  
 4. 2-(3-Phenylpropyl)pyridine