


PTC-Amino acids in Green tea

Data No. LA003-0000

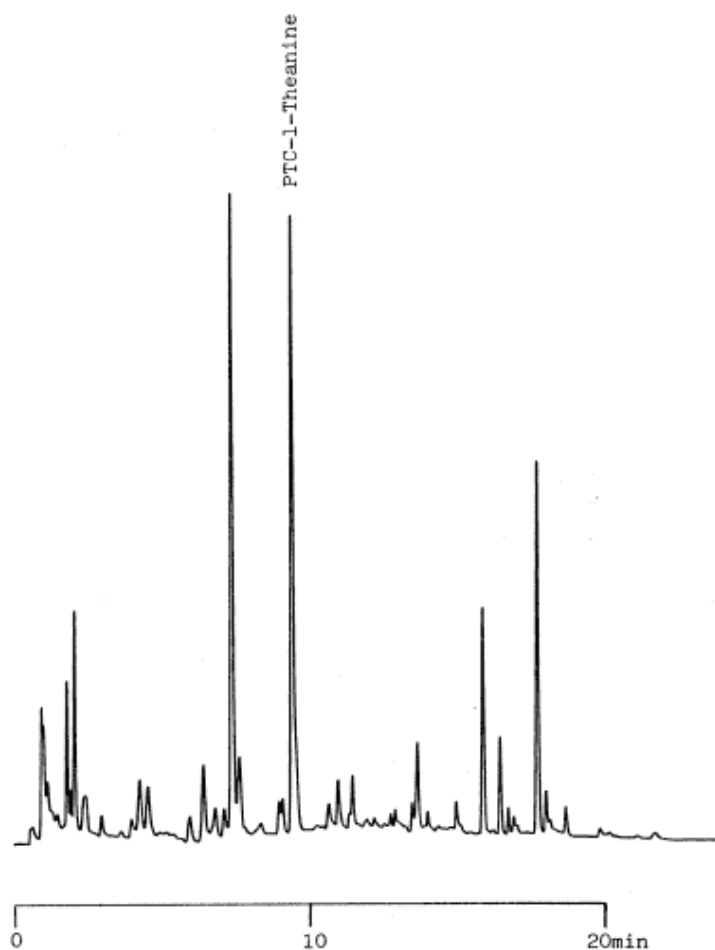

Column :Inertsil ODS-2
5 μ m
10 \times 4.0 mm I.D. (Guard) + 100 \times 4.0 mm I.D.

Column Cat. No. :5020-01101

Eluent :A: CH₃CN/0.1 M CH₃COONH₄ = 5/95
B: CH₃CN/0.1 M CH₃COONH₄ = 60/40
A/B = 100/0 - 7 min - 75/25 - 7 min - 0/100

Flow rate :1.0 mL/min.

Col. Temp. :40°C

Detection :UV 254 nm

Injection Vol. :

Sample :PTC-amino acid

Solutes :1. PTC-1-Theanine