


Free Amino acids

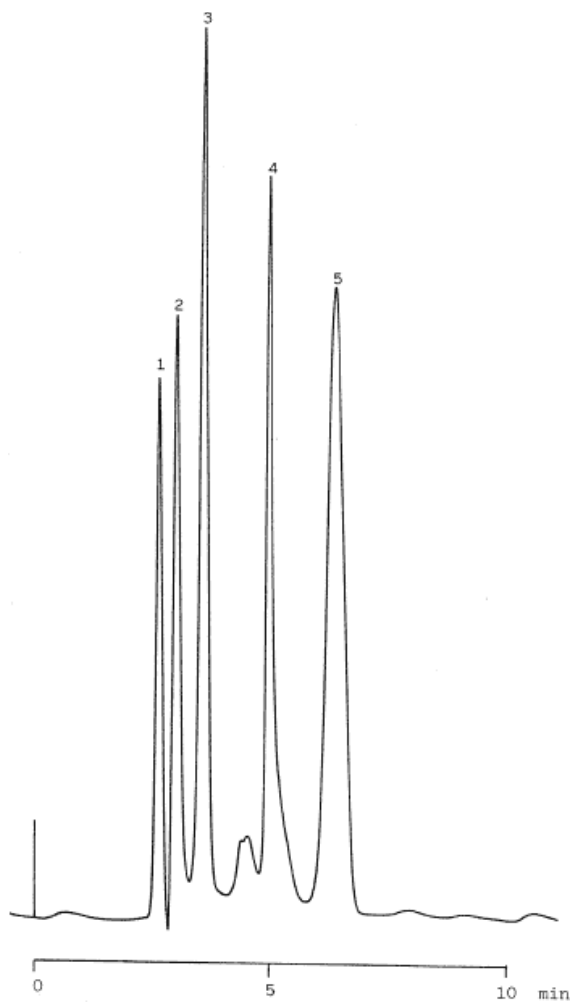

Column :Inertsil ODS
 5µm
 250 × 4.6 mm I.D.
 Column Cat. No. :5020-02102
 Eluent :20 mM KH₂PO₄ (pH 2.5; H₃PO₄)
 Flow rate :1.0 mL/min.
 Col. Temp. :40 °C
 Detection :UV 210 nm
 Injection Vol. :
 Sample :Free amino acids

Solutes :1. Lysine (0.5 mg/mL)
 2. Threonine (0.3 mg/mL)
 3. Proline (0.6 mg/mL)
 4. Valine (0.8 mg/mL)
 5. Methionine (3 µg/mL)