


InertSearch™ for LC

Inertsil® Applications

Free Amino acids

Data No. LA007-0000

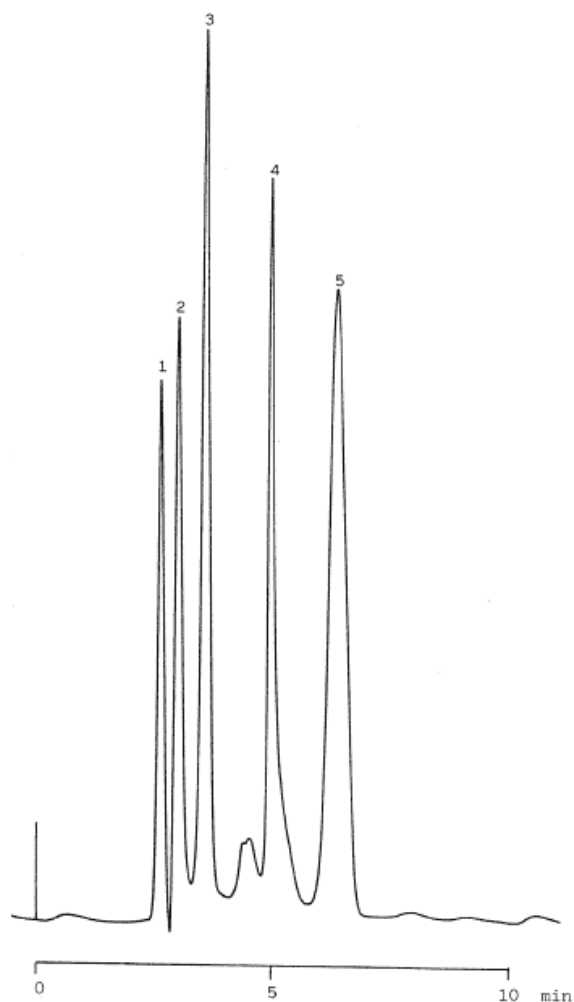

Column	:Inertsil ODS 5µm 250 × 4.6 mm I.D.	Solutes	:1. Lysine (0.5 mg/mL) 2. Threonine (0.3 mg/mL) 3. Proline (0.6 mg/mL) 4. Valine (0.8 mg/mL) 5. Methionine (3 µg/mL)
Column Cat. No.	:5020-02102		
Eluent	:20 mM KH ₂ PO ₄ (pH 2.5; H ₃ PO ₄)		
Flow rate	:1.0 mL/min.		
Col. Temp.	:40 °C		
Detection	:UV 210 nm		
Injection Vol.	:		
Sample	:Free amino acids		