


## Fat-soluble Vitamins

Data No. LA013-0000

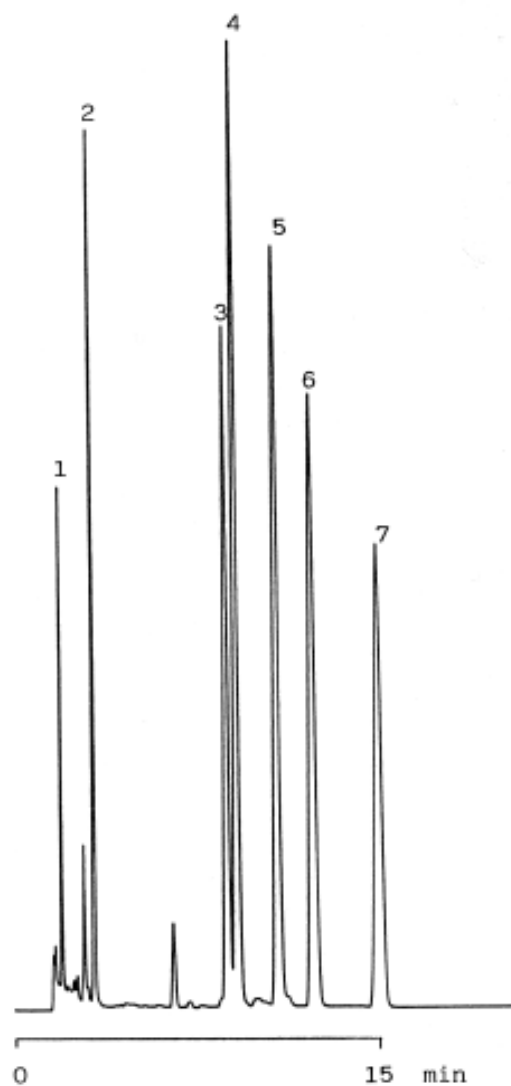

Column :Inertsil ODS-2

5µm

250 × 4.6 mm I.D

Column Cat. No. :5020-01102

Eluent :CH<sub>3</sub>CN/CH<sub>3</sub>OH = 75/25

Flow rate :1.0 mL/min.

Col. Temp. :40°C

Detection :UV 280 nm

Injection Vol. :

Sample :Fat-soluble Vitamins

Solutes

: 1. Vitamin K3

2. Vitamin A

3. Vitamin D2

4. Vitamin D3

5. Vitamin E

6. Vitamin E acetate

7. Vitamin K1