


## Vitamin E (Tocopherol) in Mother's milk

Data No. LA017-0000

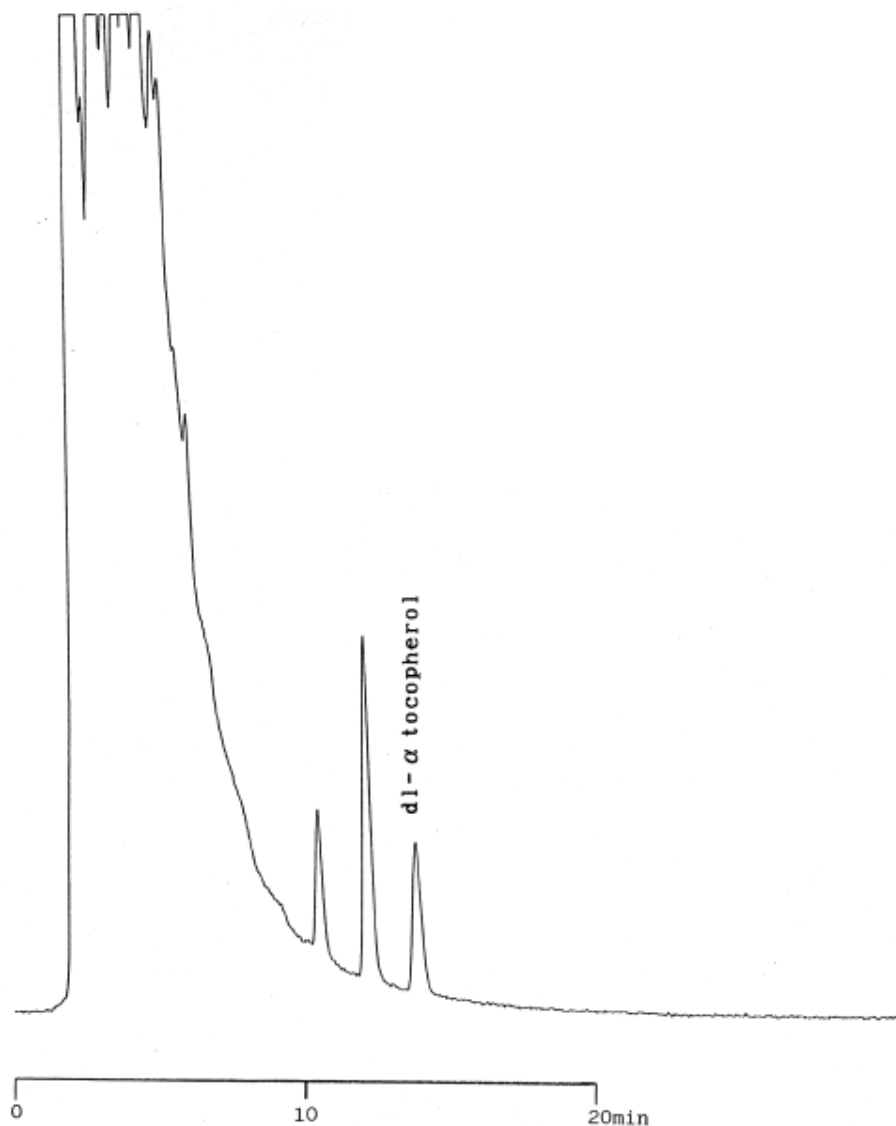

Column : Inertsil ODS-2  
5 $\mu$ m  
250  $\times$  4.6 mm I.D.

Solutes : 1.  $\alpha$ -Tocopherol

Column Cat. No.: 5020-01102

Eluent : CH<sub>3</sub>CN/CH<sub>3</sub>OH = 15/85

Flow rate : 1.5 mL/min.

Col. Temp. : Ambient

Detection : UV 295 nm

Injection Vol. :

Sample : Tocopherol