


Fat-soluble Vitamins

Data No. LA019-0000

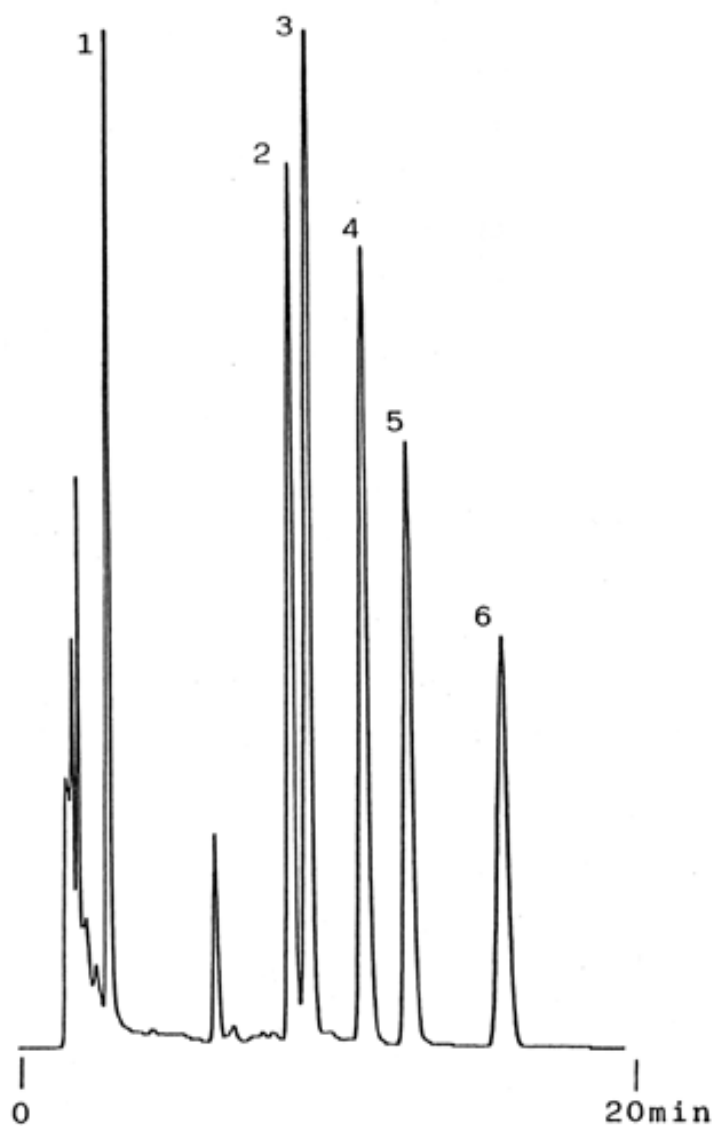

Column : Inertsil ODS-2
5µm
150 × 0.7 mm I.D.

Column Cat. No. : 5020-13073

Eluent : CH₃CN/CH₃OH = 75/25

Flow rate : 30 µL/min.

Col. Temp. :

Detection : UV 280 nm

Injection Vol. :

Sample : Fat-soluble Vitamins

Solutes : 1. Vitamin A
2. Vitamin D2
3. Vitamin D3
4. Vitamin E
5. Vitamin E acetate
6. Vitamin K1