


## Vitamins

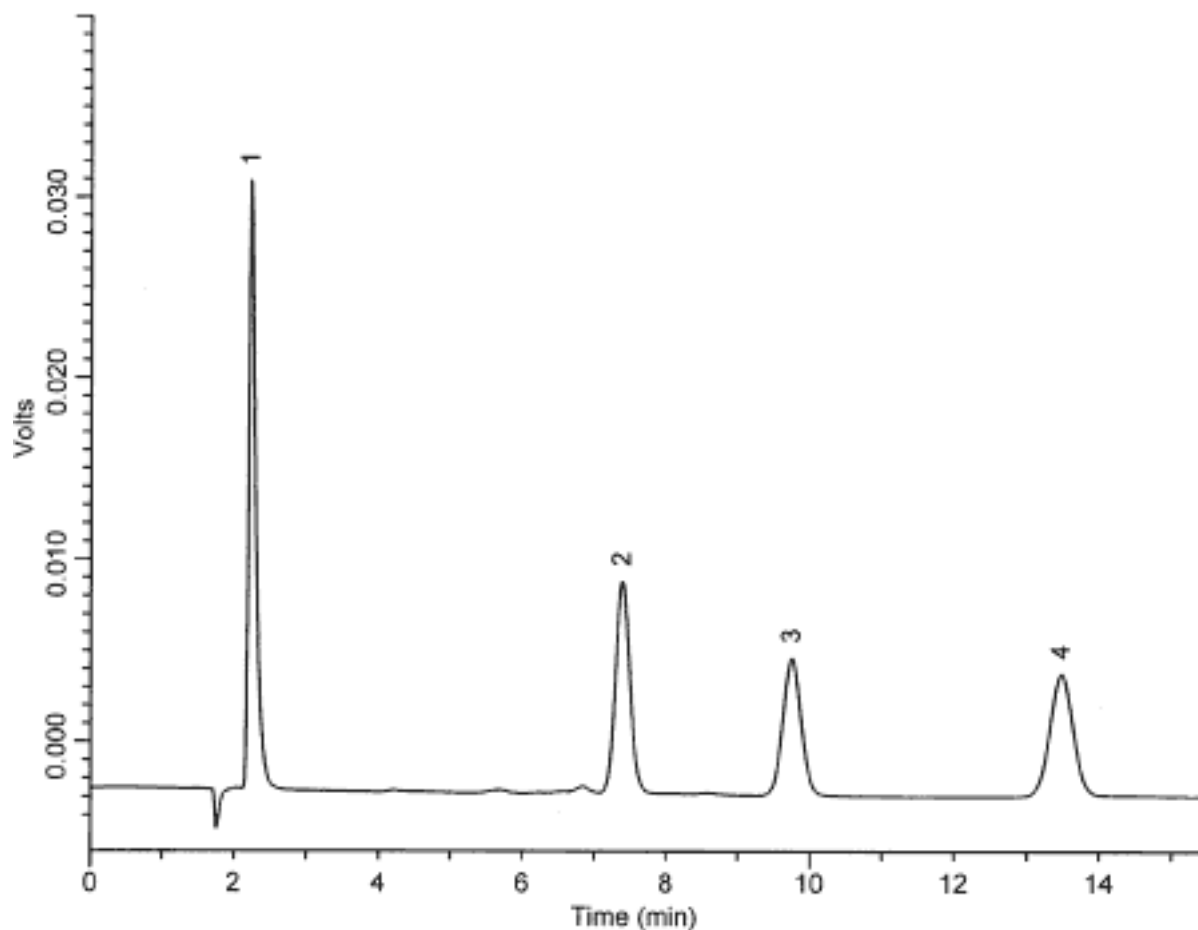

Column : Inertsil ODS-3V  
150 × 4.6 mm I.D.

Column Cat. No. :

Eluent : CH<sub>3</sub>OH

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 280 nm

Injection Vol. : 3 µL

Sample : Fat-soluble Vitamins

Solutes : 1. Vitamin K3 (0.18 mg/mL)

2. Vitamin E (0.36 mg/mL)

3. Vitamin E acetate (0.36 mg/mL)

4. Vitamin K1 (0.09 mg/mL)