


InertSearch™ for LC

Inertsil® Applications

Vitamins

Data No. LA024-0000

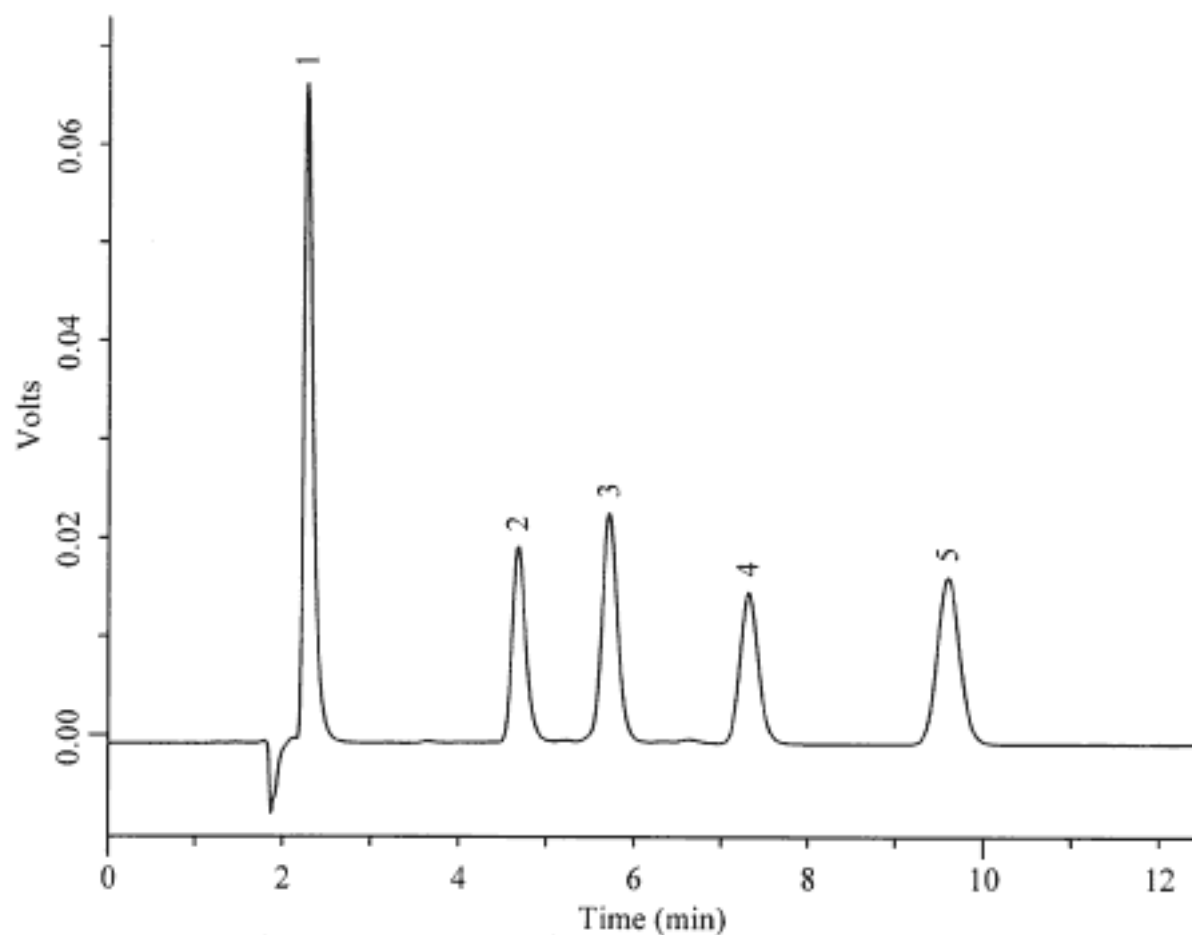

Column	: Inertsil C8-3 150 × 4.6 mm I.D.	Solutes	: 1. Vitamin K3 (0.18 mg/mL) 2. Vitamin D3 (0.03 mg/mL) 3. Vitamin E (0.36 mg/mL) 4. Vitamin E acetate (0.36 mg/mL) 5. Vitamin K1 (0.09 mg/mL)
Column Cat. No. :			
Eluent	: CH ₃ OH/H ₂ O = 95/5		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 280 nm		
Injection Vol.	: 5 µL		
Sample	: Fat-soluble Vitamins		