


# InertSearch™ for LC

Inertsil® Applications

## Drug

Data No. LA027-0000

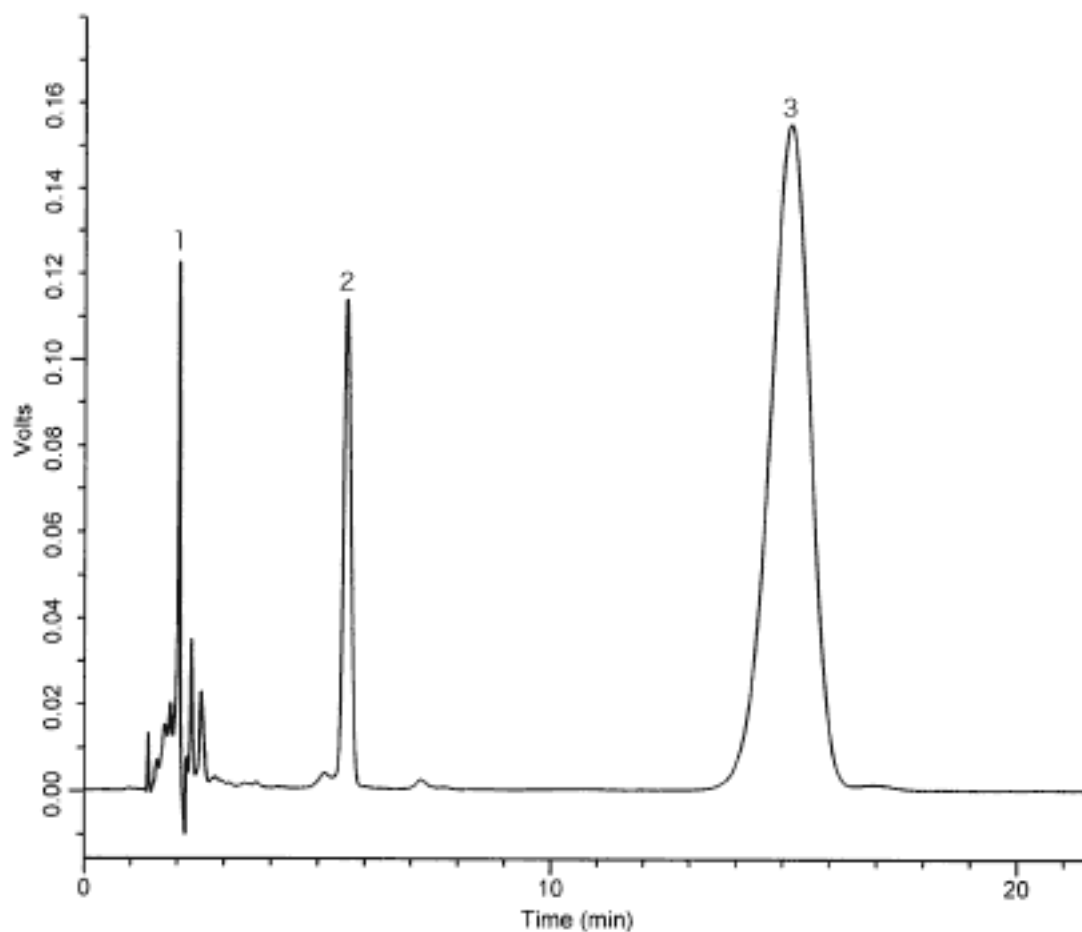

Column	: Inertsil ODS-3V 150 × 4.6 mm I.D.	Solutes	: 1. Pyridoxine (0.10 mg/mL) 2. Riboflavin (0.10 mg/mL) 3. Fursultiamine (1.00 mg/mL)
Column Cat. No. :			
Eluent	: CH <sub>3</sub> CN/0.1 % H <sub>3</sub> PO <sub>4</sub> + 5 mM IPCC-05 = 15/85		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 210 nm		
Injection Vol.	: 20 µL		
Sample	: Water-soluble Vitamins (Drugs)		