


## Vitamins

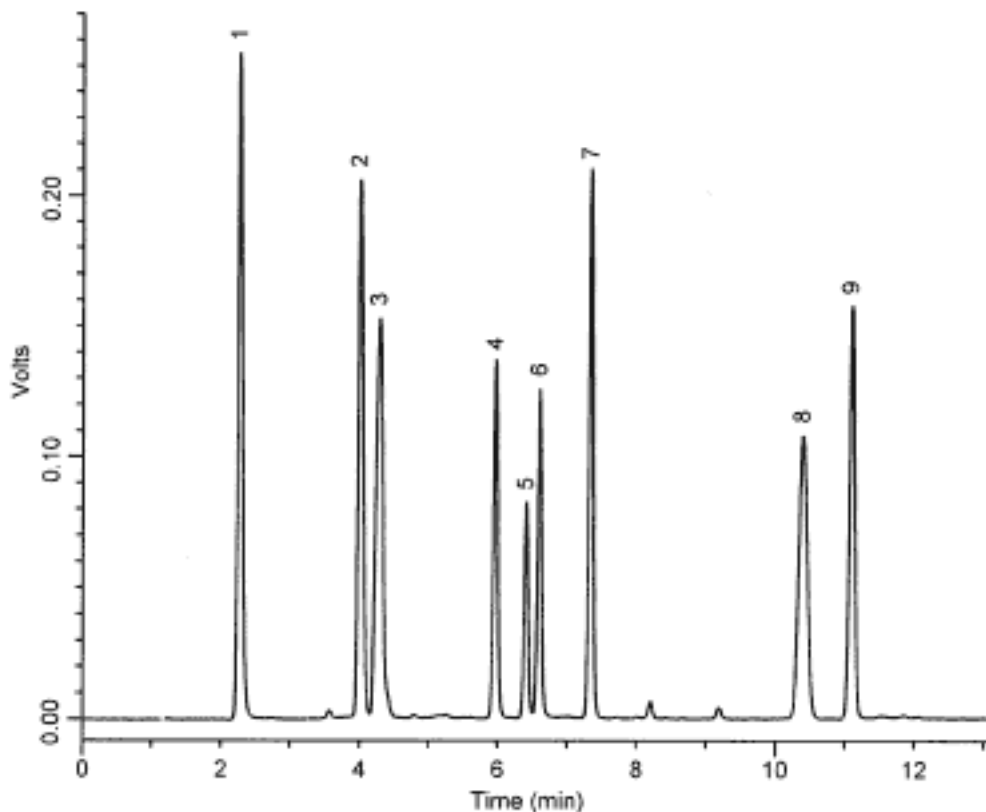

Column	: Inertsil ODS-3V 150 × 4.6 mm I.D.	Solutes	: 1. L-Ascorbic acid (0.12 mg/mL) 2. Nicotinic acid (0.10 mg/mL) 3. Nicotinamide (0.10 mg/mL) 4. Pyridoxal (0.47 mg/mL) 5. Pyridoxamine (0.37 mg/mL) 6. Pyridoxine (0.50 mg/mL) 7. Thiamine (0.07 mg/mL) 8. Cyanocobalamin (0.33 mg/mL) 9. Riboflavin (0.07 mg/mL)
Column Cat. No. :			
Eluent	: A: 0.1 % H <sub>3</sub> PO <sub>4</sub> + 5 mM IPCC-05 B: CH <sub>3</sub> CN/0.1 % H <sub>3</sub> PO <sub>4</sub> + 5 mM IPCC-05 A/B = 97/3 - 20min - 55/45		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	: 5 µL		
Sample	: Water-soluble Vitamins		