


InertSearch™ for LC

Inertsil® Applications

Vitamins

Data No. LA034-0000

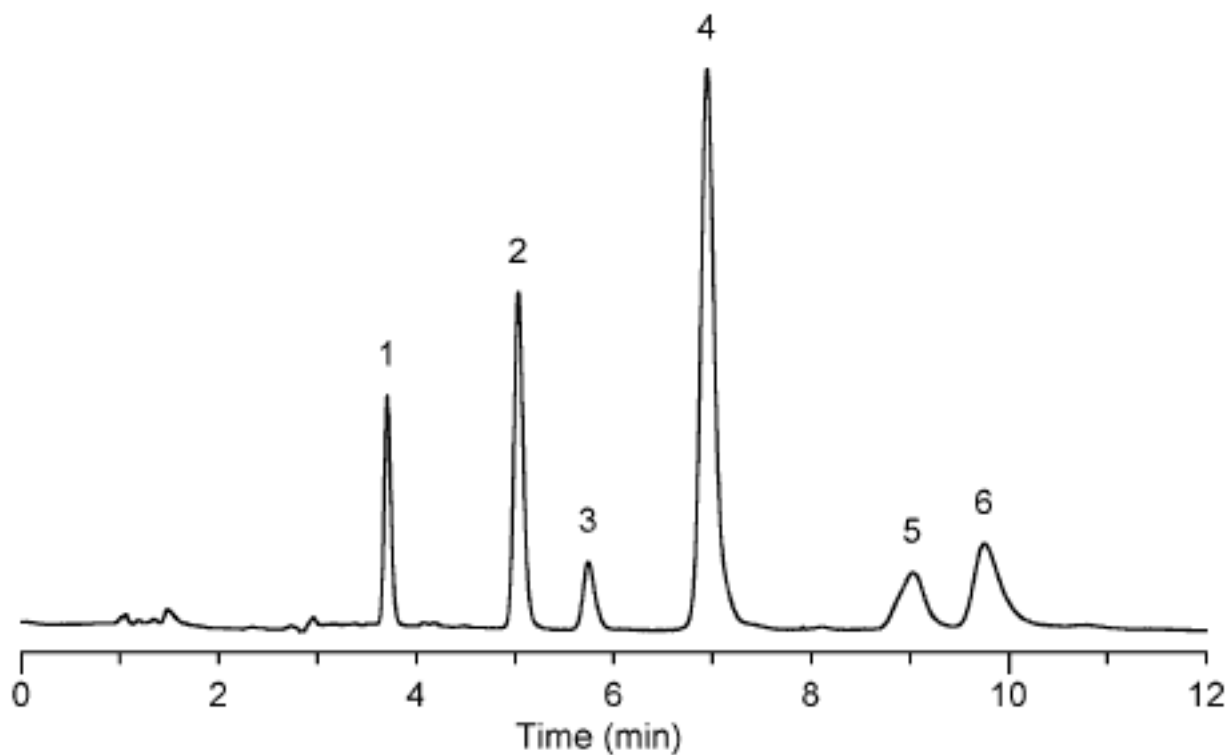

Column	: Inertsil CN-3 5µm 250 × 4.6 mm I.D.	Solutes	:1. Vitamin K3 (0.14 mg/mL, tR 3.703) 2. Vitamin A (0.14 mg/mL, tR 5.033) 3. Vitamin A acetate (0.14 mg/mL, tR 5.740) 4. Vitamin D3 (0.14 mg/mL, tR 6.943) 5. Vitamin E (0.14 mg/mL, tR 9.027) 6. Vitamin K1 (0.14 mg/mL, tR 9.750)
Column Cat. No.:			
Eluent	: CH ₃ CN/H ₂ O = 60/40		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 280 nm		
Injection Vol.	: 1 µL		
Sample	: Fat-soluble Vitamins		