


InertSearch™ for LC

Inertsil® Applications

Vitamins

Data No. LA035-0000

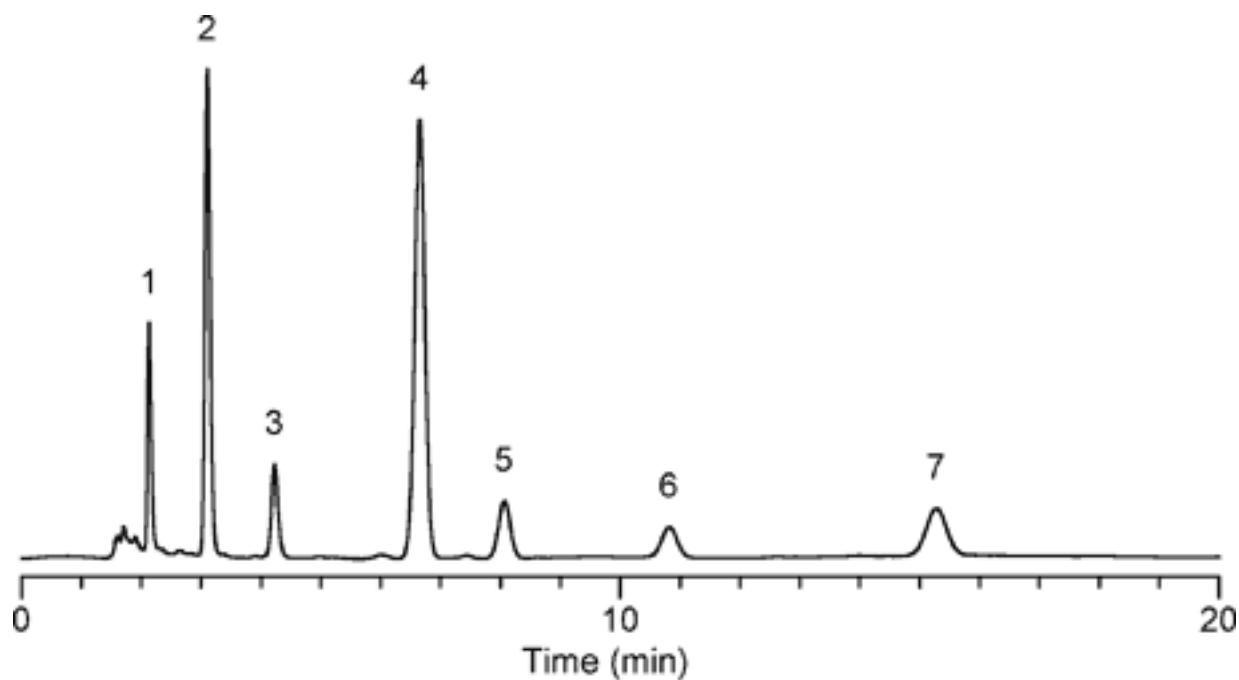

Column	: Inertsil ODS-3V 5 μ m 150 \times 4.6 mm I.D.	Solutes	: 1. Vitamin K3 (0.14 mg/mL, tR 2.130) 2. Vitamin A (0.14 mg/mL, tR 3.103) 3. Vitamin A acetate (0.14 mg/mL, tR 4.227) 4. Vitamin D3 (0.14 mg/mL, tR 6.650) 5. Vitamin E (0.14 mg/mL, tR 8.063) 6. Vitamin E acetate (0.14 mg/mL, tR 10.820) 7. Vitamin K1 (0.14 mg/mL, tR 15.270)
Column Cat. No.	:		
Eluent	: CH ₃ OH		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 280 nm		
Injection Vol.	: 1 μ L		
Sample	: Fat-soluble Vitamins		