


Vitamins

Data No. LA038-0000

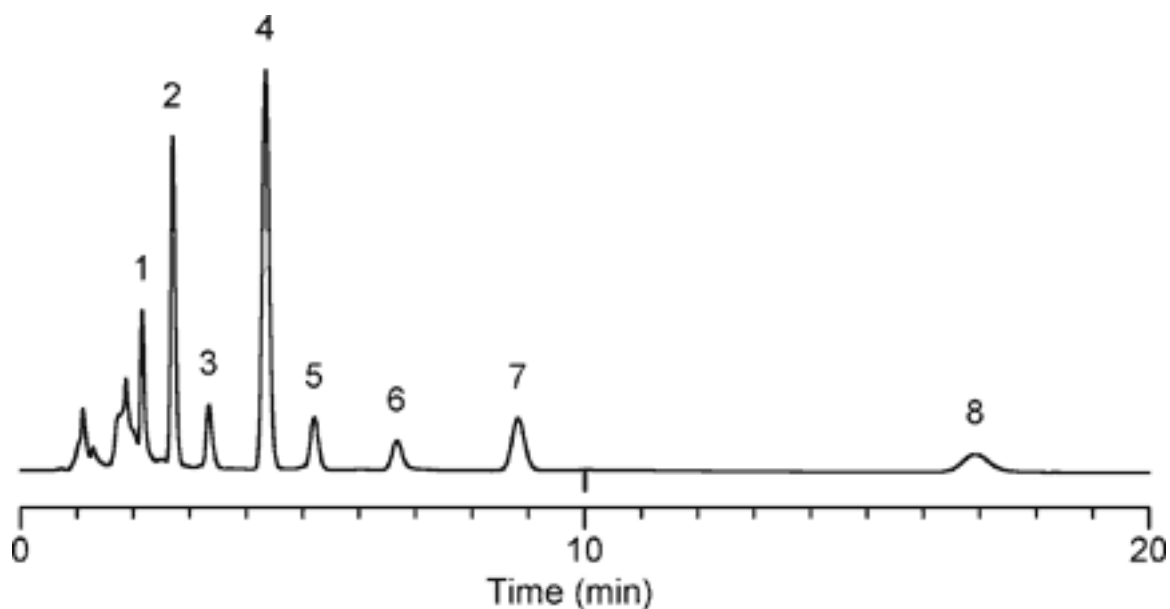

Column	: Inertsil C8-3 5µm 150 × 4.6 mm I.D.	Solutes	: 1. Vitamin K3 (0.07 mg/mL, tR 2.153) 2. Vitamin A (0.07 mg/mL, tR 2.697) 3. Vitamin A acetate (0.07 mg/mL, tR 3.333) 4. Vitamin D3 (0.07 mg/mL, tR 4.343) 5. Vitamin E (0.07 mg/mL, tR 5.203) 6. Vitamin E acetate (0.07 mg/mL, tR 6.670) 7. Vitamin K1 (0.07 mg/mL, tR 8.817) 8. β-Carotene (0.13 mg/mL, tR 16.920)
Column Cat. No. :			
Eluent	: CH ₃ OH/H ₂ O = 95/5		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: 0 - 10 min: UV 280 nm, 10 - 20 min: VIS 450 nm		
Injection Vol.	: 3 µL		
Sample	: Fat-soluble Vitamins		