


Analysis of Water-soluble Vitamins

Data No. LA047-0000

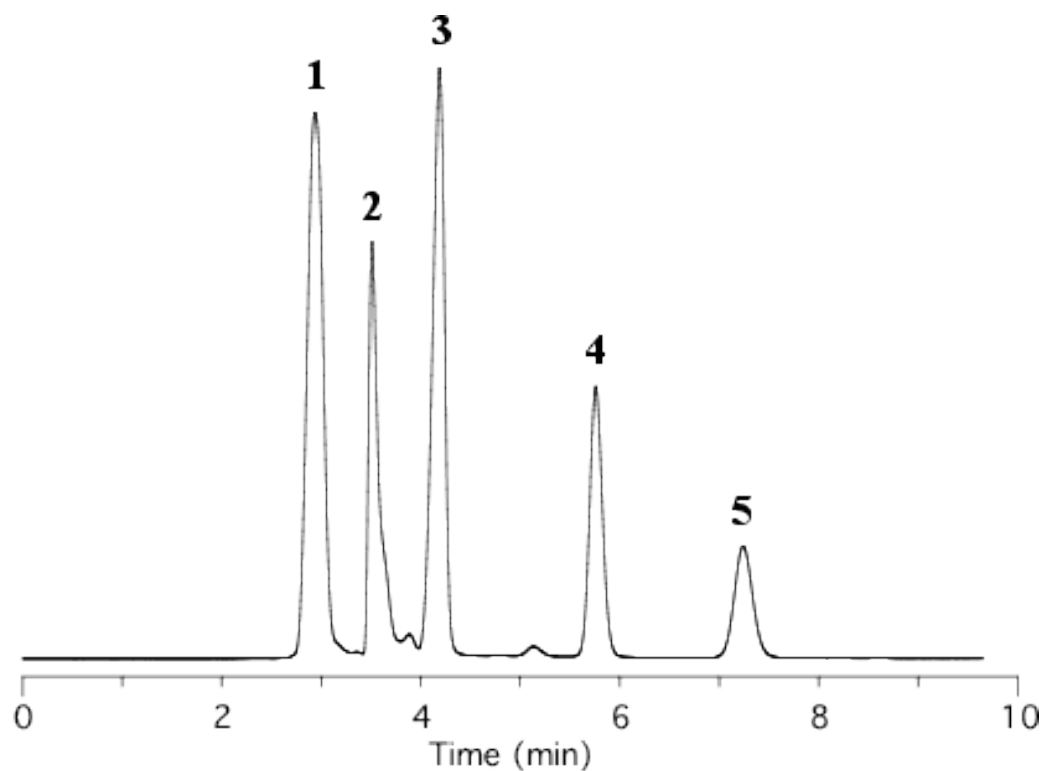

Column	: Inertsil Diol 5 μ m 250 x 4.6 mm I.D.	Solutes	: 1. Vitamin B1 (Thiamine) 2. Vitamin B3 (Niacin, Nicotinic acid) 3. Vitamin B3 (Nicotinamide) 4. Vitamin C (L-Ascorbic acid) 5. Vitamin B2 (Riboflavin)
Column Cat. No. :			
Eluent	: A) CH ₃ CN B) H ₂ O C) TFA A/B/C = 90/10/0.1, w/w/v		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		