


## Analysis of Carotenoids

Data No. LA049-0000

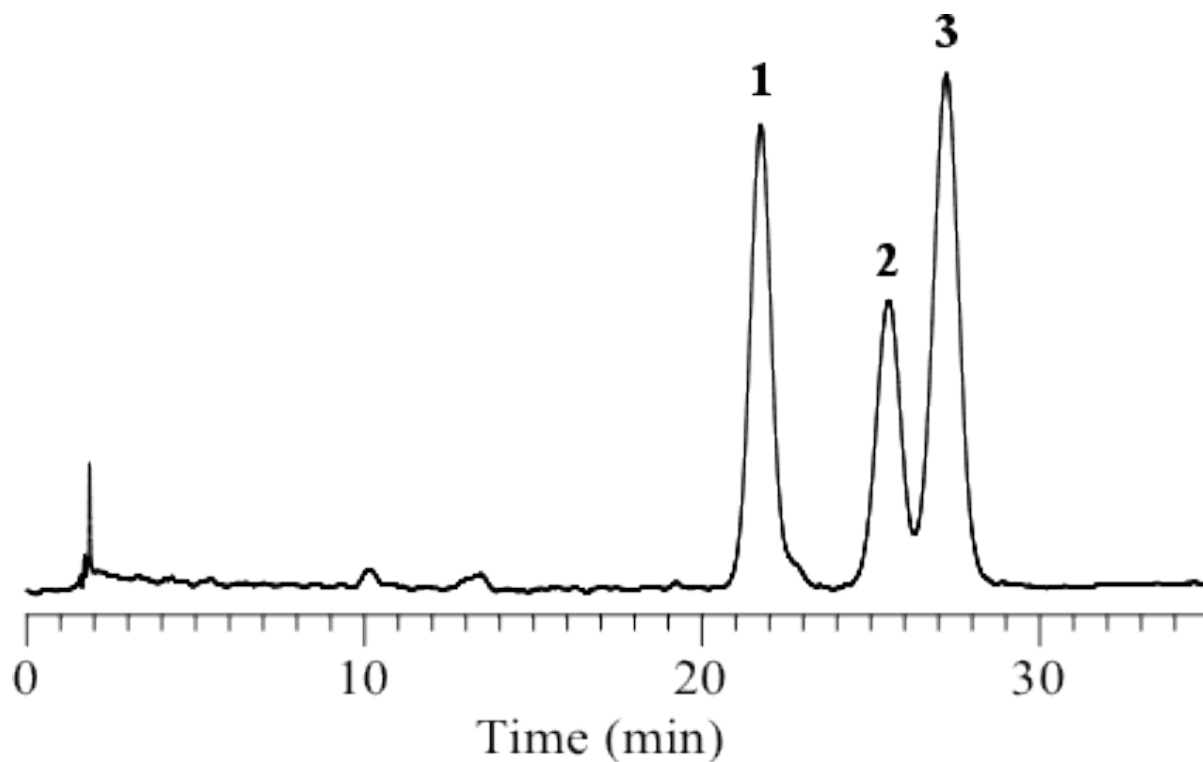

Column	: Inertsil ODS-3	Solutes	: 1. Lycopene
	5 $\mu$ m		2. $\alpha$ -Carotene
	150 x 4.6 mm I.D.		3. $\beta$ -Carotene

Column Cat. No.: 5020-01731

Eluent : A) CH<sub>3</sub>OH  
B) THF  
A/B = 95/5, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : VIS 450 nm

Injection Vol. :

Sample :