


# InertSearch™ for LC

Inertsil® Applications

## Analysis of Carotenoids

Data No. LA050-0000

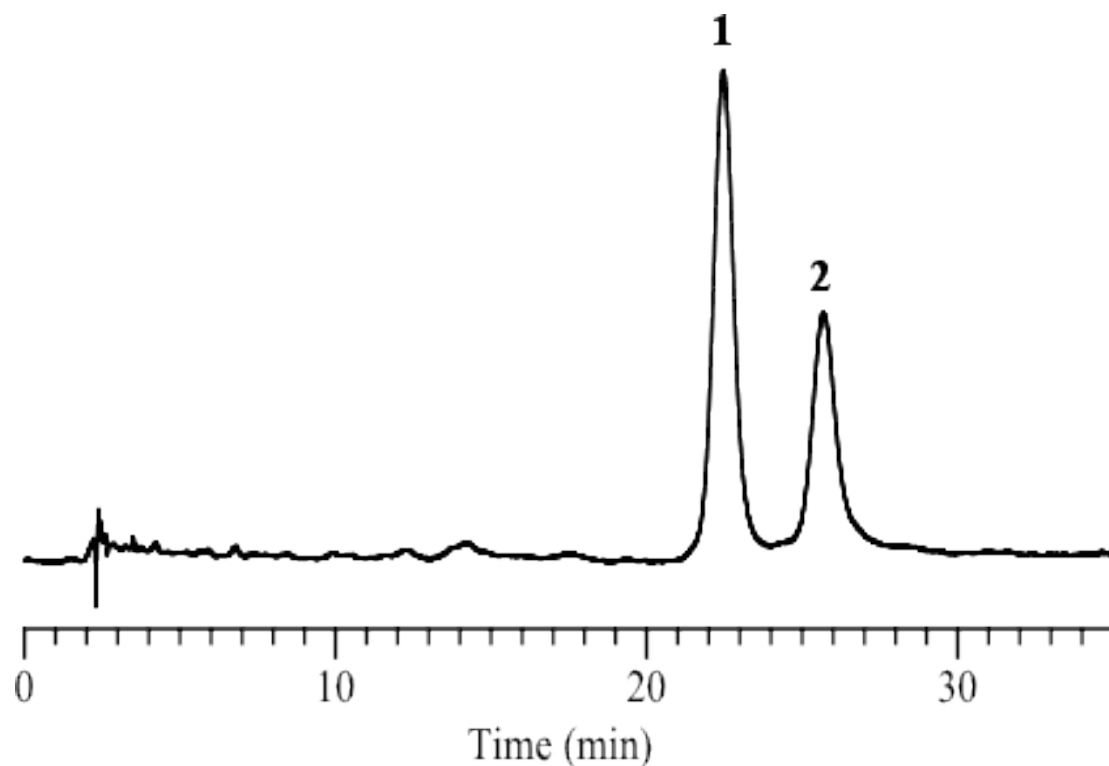

Column : Inertsil ODS-P  
5  $\mu$ m  
250 x 4.6 mm I.D.

Solutes : 1.  $\alpha$ -Carotene  
2.  $\beta$ -Carotene

Column Cat. No. : 5020-02002

Eluent : A) CH<sub>3</sub>OH  
B) THF  
A/B = 95/5, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 30°C

Detection : VIS 450 nm

Injection Vol. :

Sample :