


1. Food additives

Data No. LA056-0000

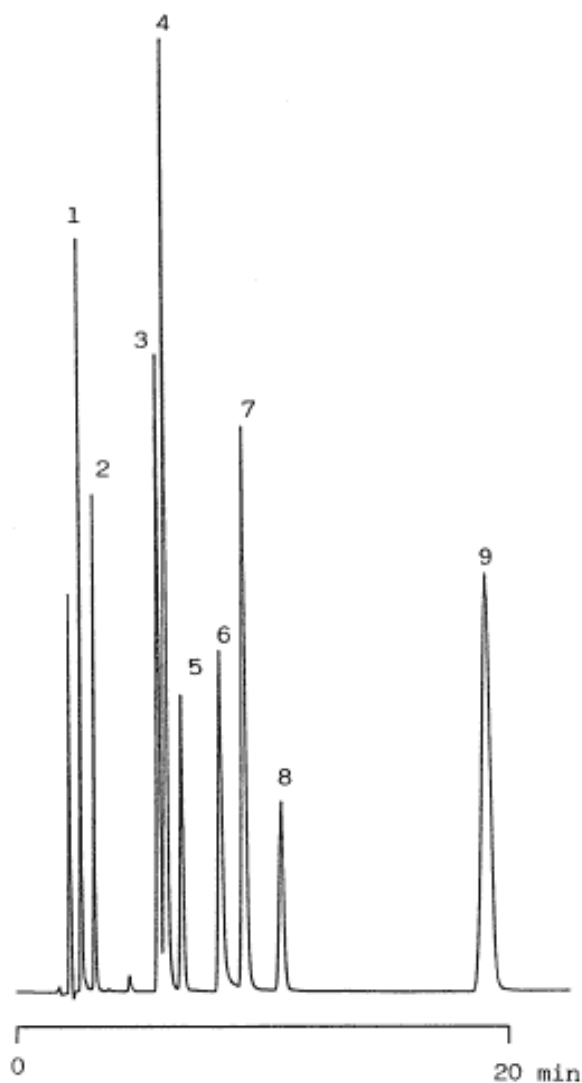

Column : Inertsil Ph
5µm
150 × 4.6 mm I.D.

Column Cat. No. : 5020-01301

Eluent : CH₃CN/50 mM KH₂PO₃ + 0.1 % H₃PO₄ = 25/75

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. :

Sample : Food Additives

Solutes : 1. Sodium saccharin
2. p-Hydroxy benzoic acid
3. Sorbic acid
4. Benzoic acid
5. p-Hydroxy benzoic acid methyl ester
6. Dehydroacetic acid
7. p-Toluic acid
8. p-Hydroxy benzoic acid ethyl ester
9. p-Hydroxy benzoic acid n-propyl ester