


5. Ginsenosides

Data No. LA060-0000

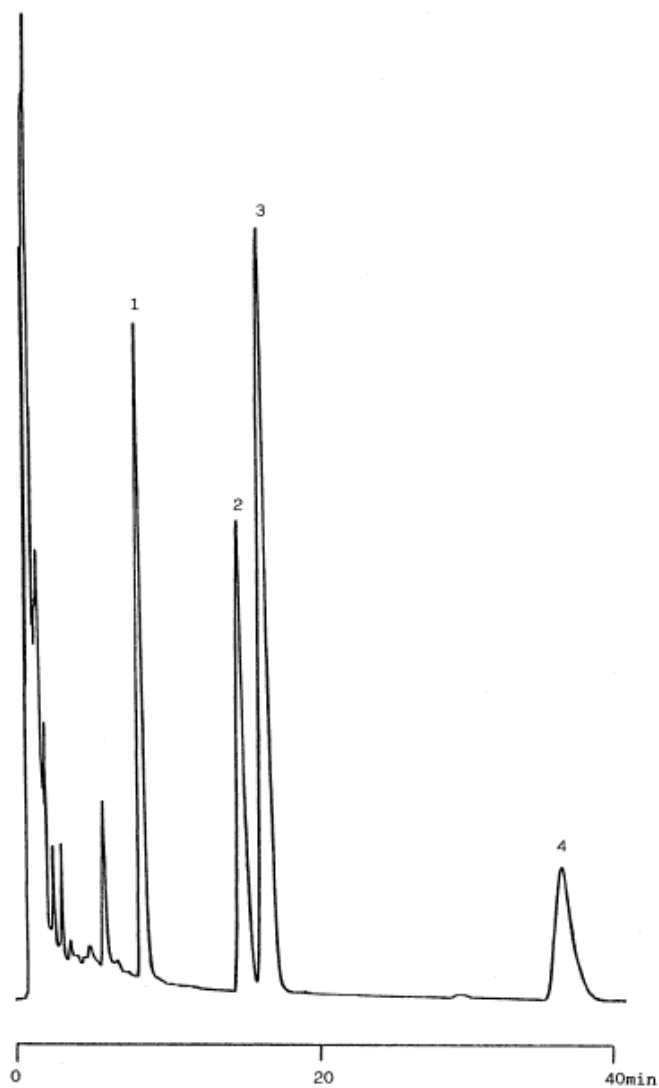

Column : Inertsil ODS
5 μ m
150 \times 4.6 mm I.D.

Column Cat. No. : 5020-02101

Eluent : CH₃CN/H₂O = 27/73

Flow rate : 1.0 mL/min.

Col. Temp. : Ambient

Detection : UV 203 nm

Injection Vol. :

Sample : Ginsenosides

Solutes : 1. Ginsenoside Rb1
2. Ginsenoside Rc
3. Ginsenoside Rb2
4. Ginsenoside Rd