


6. Synthetic food dyes

Data No. LA061-0000

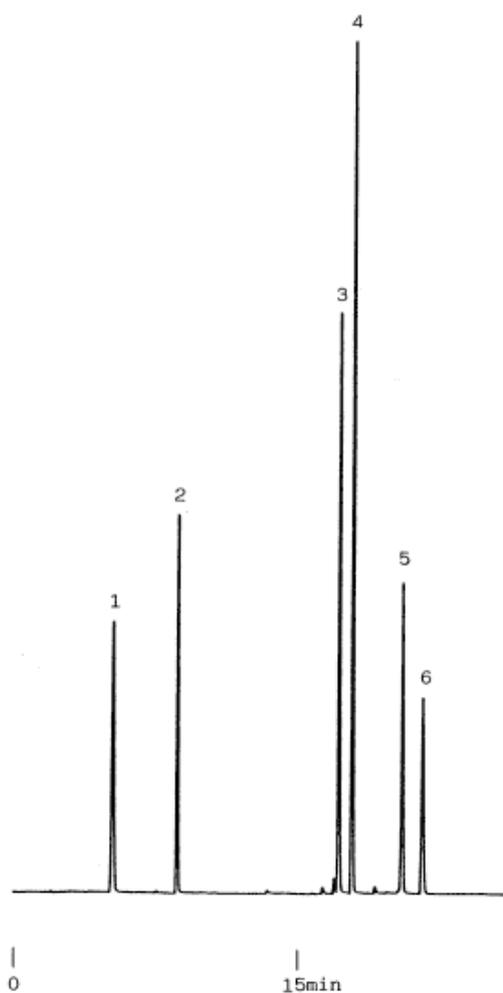

Column :Inertsil ODS-2
5µm
150 × 4.6 mm I.D.

Column Cat. No. :5020-01101

Eluent :A: CH₃OH/10 mM (NH₄)₂HPO₄ = 10/90
B: CH₃OH/10 mM (NH₄)₂HPO₄ = 70/30
A/B = 100/0 - 15 min - 0/100

Flow rate :1.0 mL/min.

Col. Temp. :50°C

Detection :VIS 530 nm

Injection Vol. :

Sample :Synthetic Food Dyes

Solutes : 1. R-2
2. R-102
3. R-3
4. R-106
5. R-104
6. R-105
each 100 ng