


11. TFA-Aflatoxins

Data No. LA066-0000

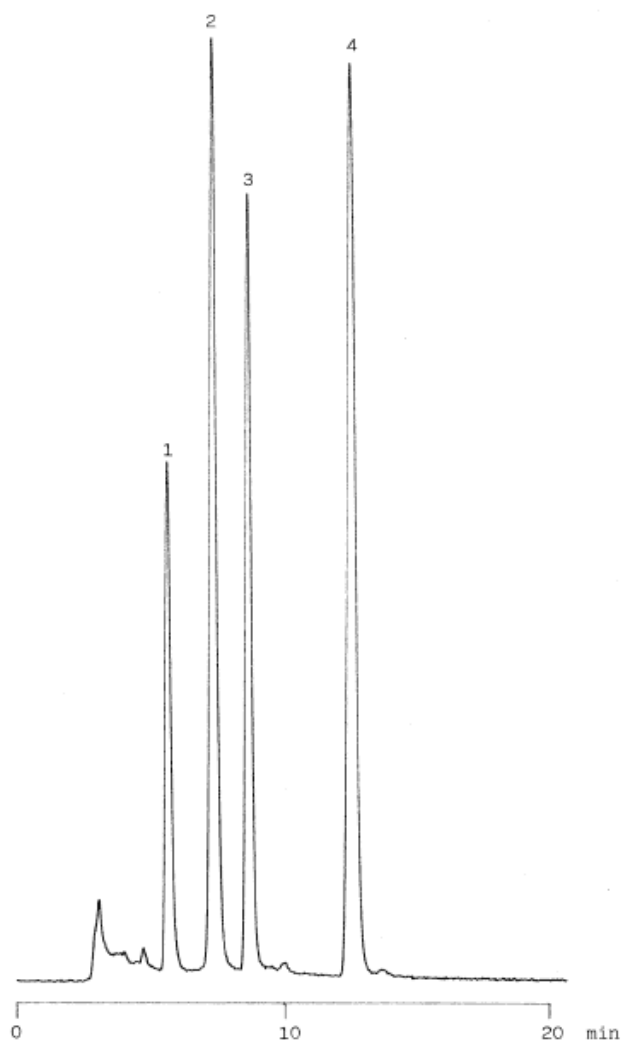

Column : Inertsil ODS-2
5 μ m
150 \times 4.6 mm I.D.

Column Cat. No. : 5020-01101

Eluent : CH₃OH/H₂O = 45/55

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : FL Ex 360 nm
Em 450 nm

Injection Vol. :

Sample : TFA-Aflatoxins

Solutes : 1. TFA-Aflatoxin G1
2. TFA-Aflatoxin B1
3. TFA-Aflatoxin G2
4. TFA-Aflatoxin B2