


## 22. Analysis of alpha-Trenbolone

Data No. LA077-0000

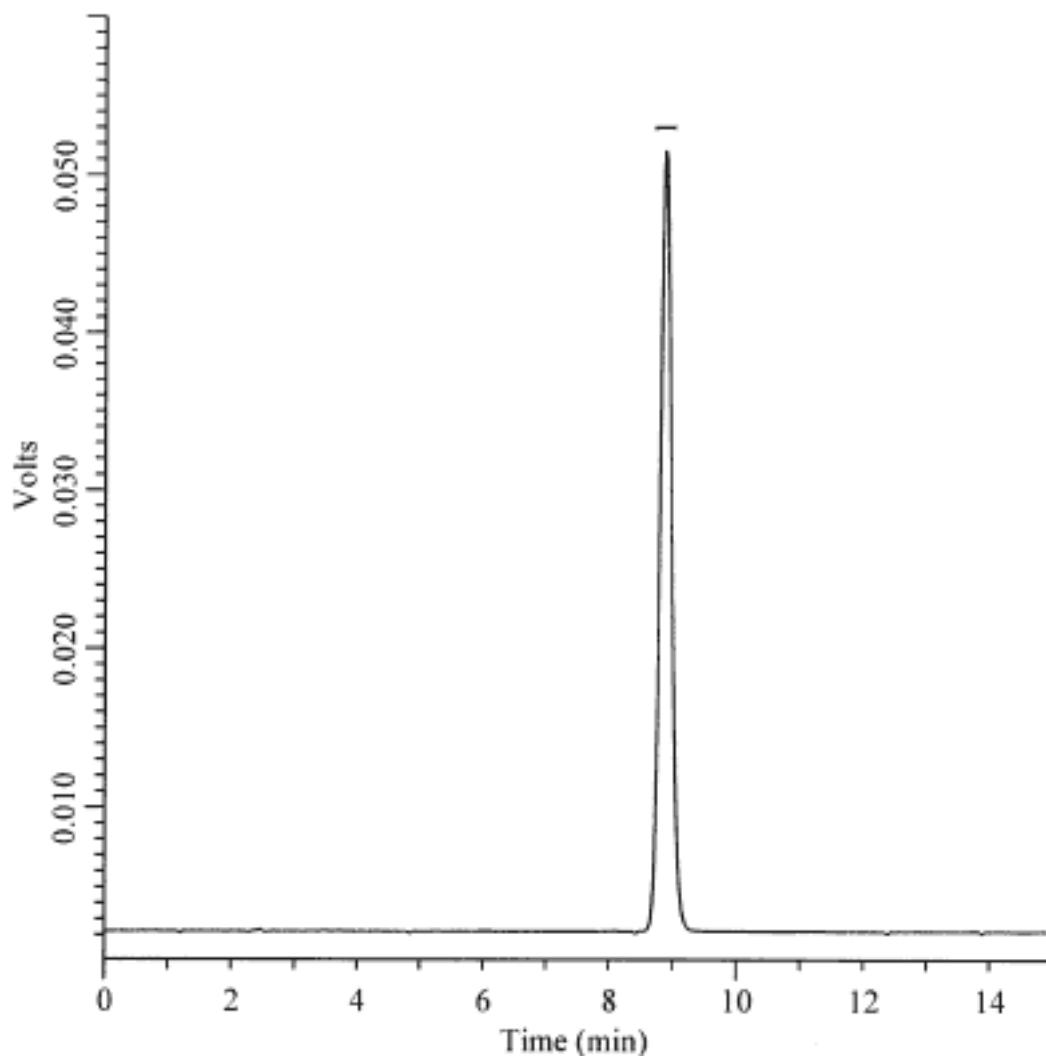

Column : Inertsil ODS-3V  
150 × 4.6 mm I.D.

Column Cat. No. :

Eluent : CH<sub>3</sub>CN/H<sub>2</sub>O = 5/6

Flow rate : 0.6 mL/min.

Col. Temp. : 50°C

Detection : UV 340 nm

Injection Vol. : 1 µL

Sample : α-Trenbolone

Solutes : 1. α-Trenbolone (100 µg/mL)