


## 23. Analysis of beta-Trenbolone

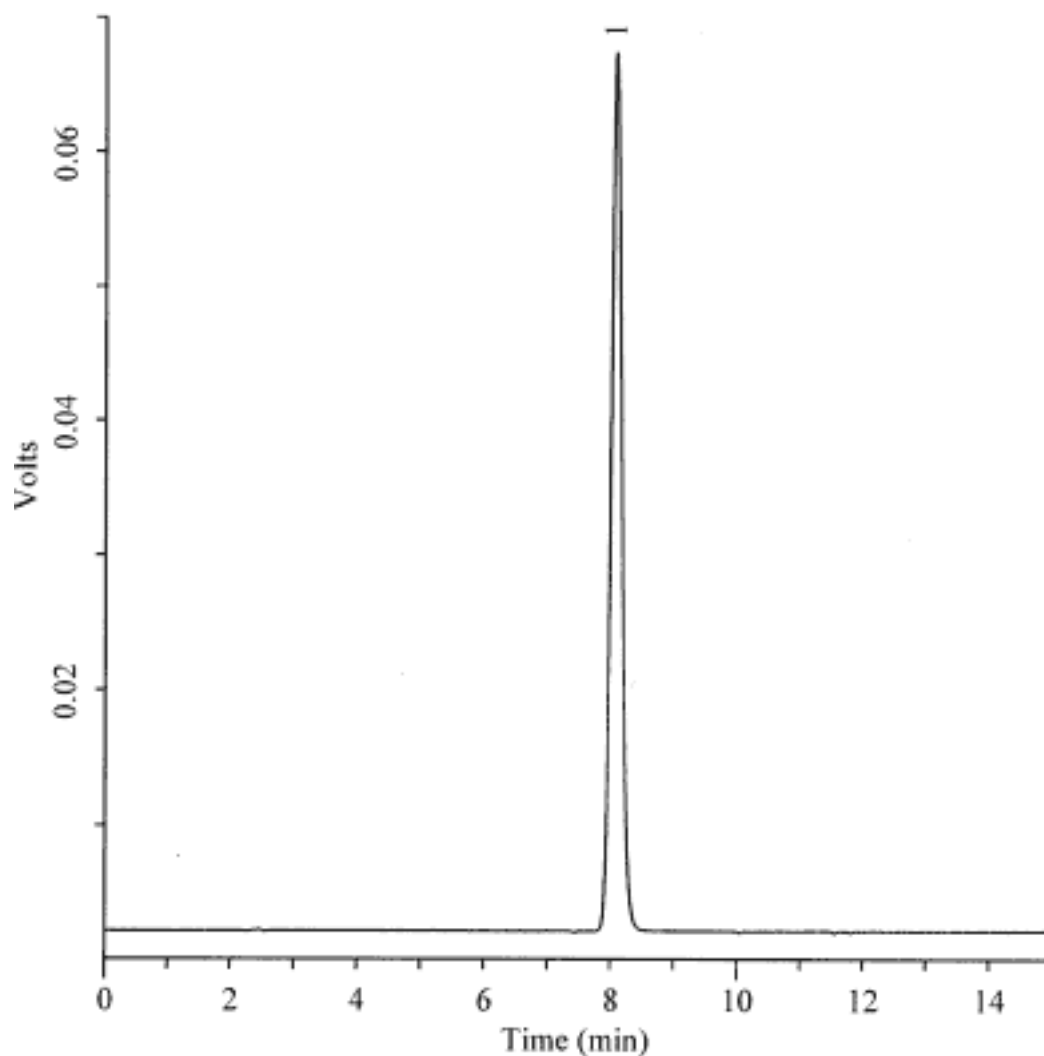

Column : Inertsil ODS-3V  
           150 × 4.6 mm I.D.  
 Column Cat. No. :  
 Eluent : CH<sub>3</sub>CN/H<sub>2</sub>O = 5/6  
 Flow rate : 0.6 mL/min.  
 Col. Temp. : 50°C  
 Detection : UV 340 nm  
 Injection Vol. : 1 µL  
 Sample : β - Trenbolone  
 Solutes : 1. β - Trenbolone (100 µg/mL)