


28. Analysis of Tetracyclines

Data No. LA083-0000

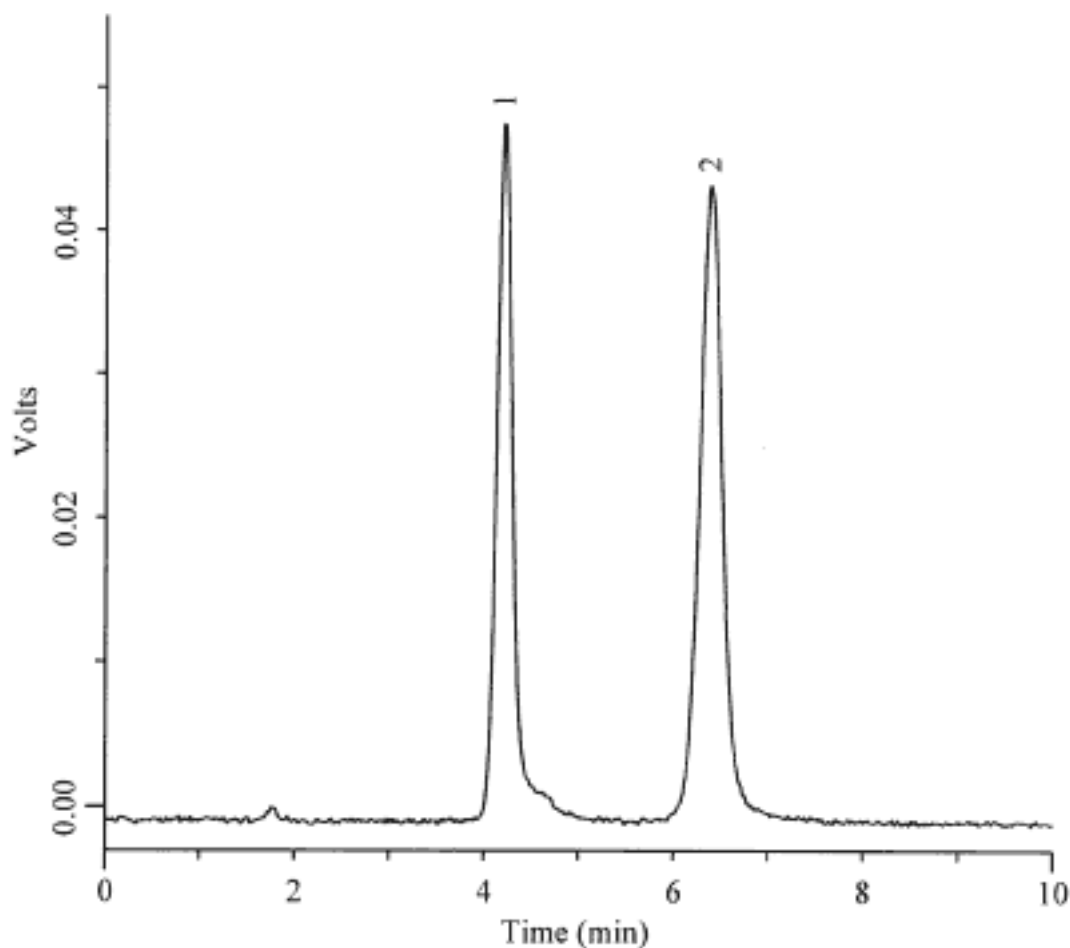

Column : Inertsil ODS-3V
150 × 4.6 mm I.D.

Solutes : 1. Oxytetracycline (100 µg/mL)
2. Tetracycline (100 µg/mL)

Column Cat. No. :

Eluent : Imidazole buffer (pH 7.2)/CH₃OH = 80/20

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : FL Ex 380 nm
Em 520 nm

Injection Vol. : 1 µL

Sample : Tetracyclines