


29. Analysis of Hormones

Data No. LA084-0000

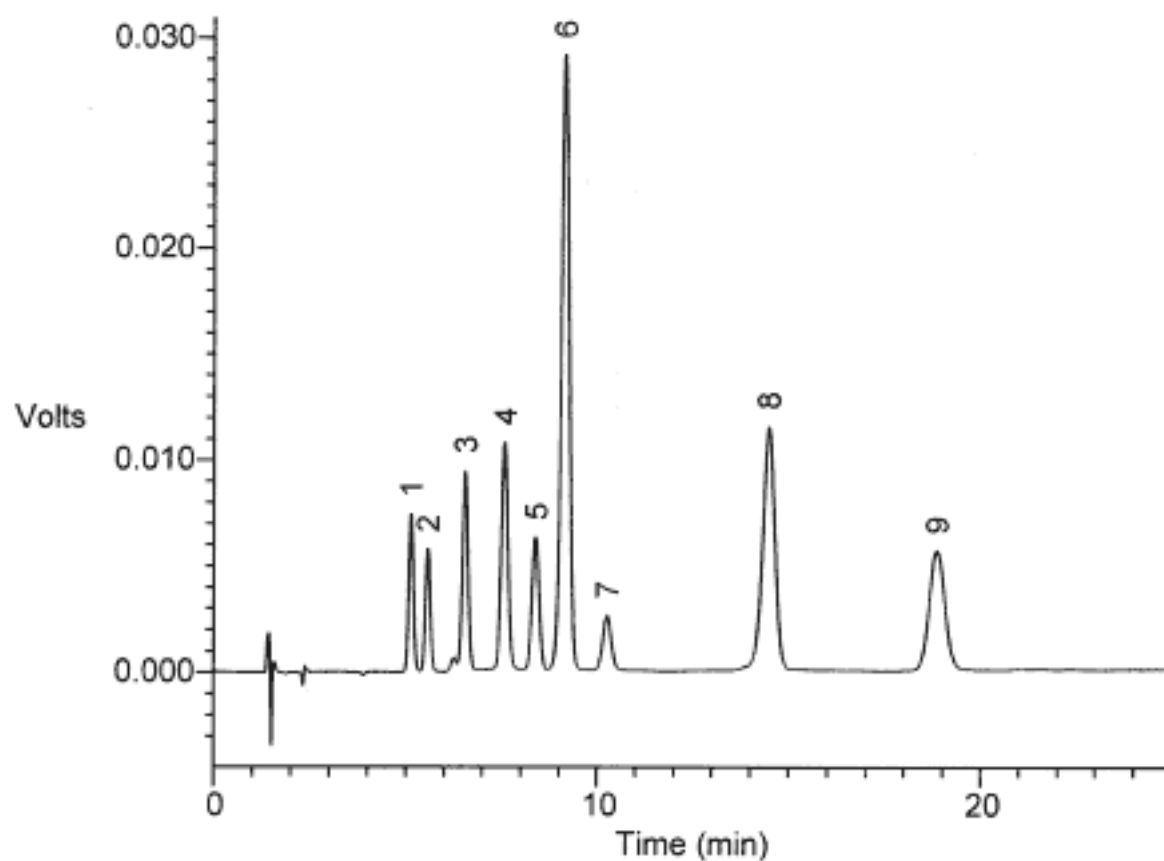

Column : Inertsil ODS-3V
150 × 4.6 mm I.D.

Column Cat. No. :

Eluent : CH₃CN/THF/H₂O = 35/8/57

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 230 nm

Injection Vol. : 5 µL

Sample : Hormones

Solutes : 1. β-Trenbolone (50 mg/L)
2. α-Trenbolone (50mg/L)
3. β-Zearalenol (30 mg/L)
4. Testosterone (50 mg/L)
5. Zeranol (50 mg/L)
6. α-Zearalenol (30 mg/L)
7. α-Estradiol (50 mg/L)
8. Zearalenone (50 mg/L)
9. Diethylstilbestrol (50 mg/L)