


## 31. Analysis of Aflatoxin

Data No. LA086-0000

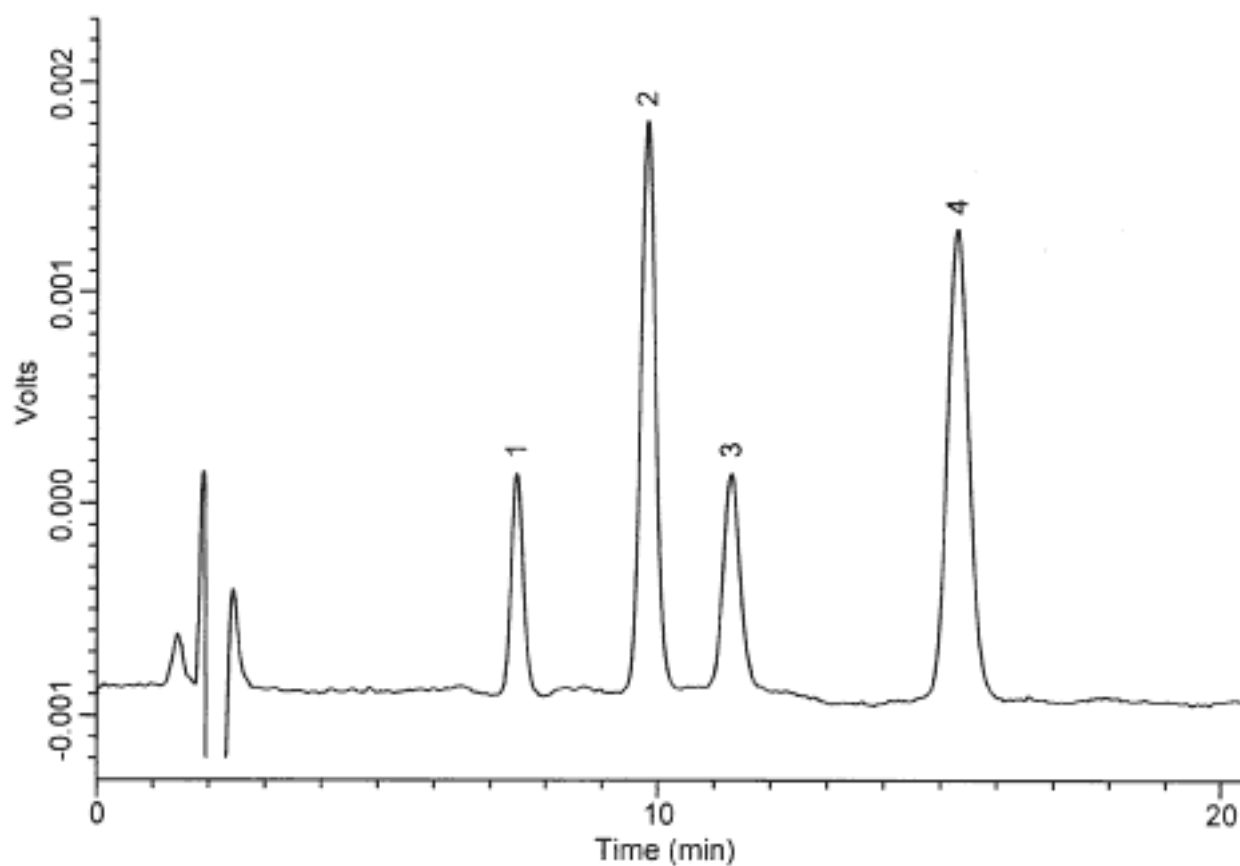

Column	:Inertsil C8-3 150 × 4.6 mm I.D.	Solutes	: 1. Aflatoxin G2 (0.06 µg/mL) 2. Aflatoxin G1 (0.20 µg/mL) 3. Aflatoxin B2 (0.06 µg/mL) 4. Aflatoxin B1 (0.20 µg/mL)
Column Cat. No. :			
Eluent	:CH <sub>3</sub> OH/H <sub>2</sub> O = 40/60		
Flow rate	:1.0 mL/min.		
Col. Temp.	:40°C		
Detection	:UV 365 nm		
Injection Vol.	:100 µL		
Sample	:Aflatoxins		