


41. Sugar alcohol

Data No. LA095-0000

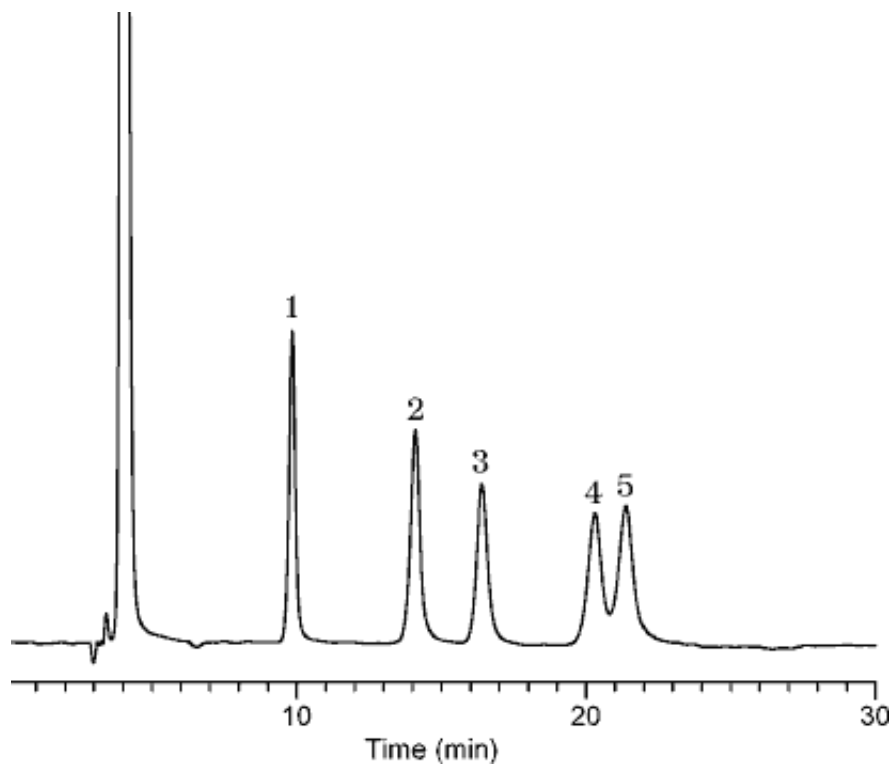

| | | | |
|-----------------|-----------------------------------------------|---------|------------------------------------------------------------------------------------------------------------------------------|
| Column | : Inertsil NH2 5µm 250 × 4.6 mm I.D. | Solutes | : 1. Erythritol (2 mg/mL) 2. Xylitol (2 mg/mL) 3. Fructose (6 mg/mL) 4. Sorbitol (2 mg/mL) 5. Mannitol (2 mg/mL) |
| Column Cat. No. | : 5020-05546 | | |
| Eluent | : CH ₃ CN/H ₂ O = 85/15 | | |
| Flow rate | : 1.0 mL/min. | | |
| Col. Temp. | : 30°C | | |
| Detection | : RI | | |
| Injection Vol. | : | | |
| Sample | : Sugar Alcohols | | |