


InertSearch™ for LC

Inertsil® Applications

43. Sugars

Data No. LA097-0000

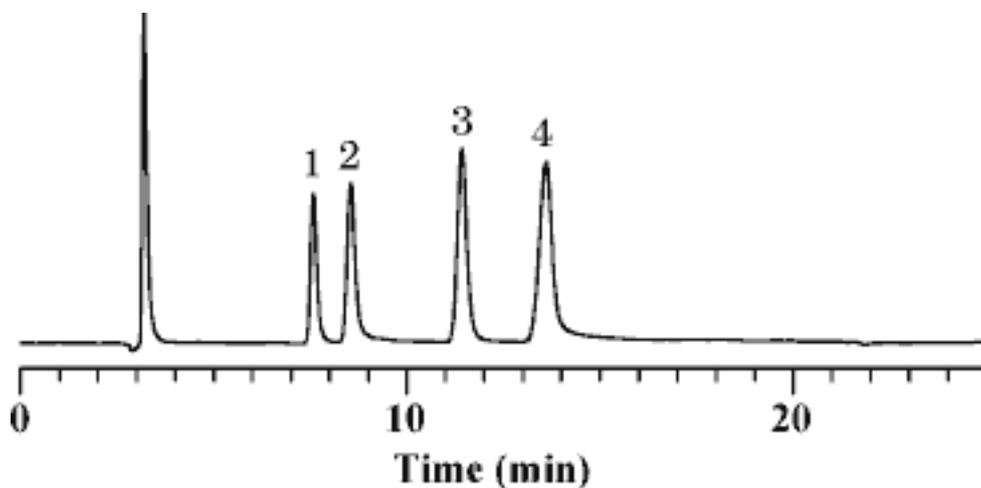

Column : Inertsil NH2
5 μ m
250 \times 4.6 mm I.D.

Column Cat. No. : 5020-05546

Eluent : CH₃CN/H₂O = 75/25

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : RI

Injection Vol. : 3 μ L

Sample : Sugars

Solutes : 1. Fructose (7.5 mg/mL)
2. Glucose (7.5 mg/mL)
3. Sucrose (7.5 mg/mL)
4. Maltose (7.5 mg/mL)