


Analysis of Parabens

Data No. LA103-0000

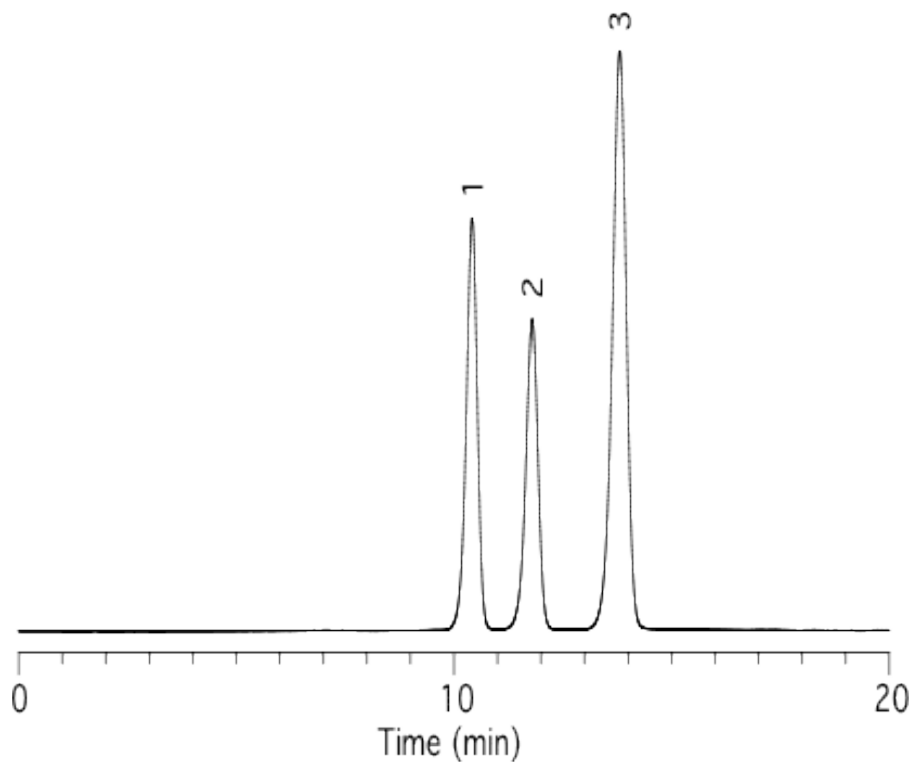

Column	: Inertsil Diol	Solutes	: 1. n-Butyl paraben
	5 µm		2. n-Propyl paraben
	150 x 4.6 mm I.D.		3. Ethyl paraben

Column Cat. No. : 5020-05645

Eluent : A) n-Hexane
B) Ethanol
A/B = 95/5, w/w

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :